

# National Human Exposure Assessment Survey (NHEXAS)

## *Maryland Study*

## Quality Systems and Implementation Plan for Human Exposure Assessment

Emory University  
Atlanta, GA 30322

Cooperative Agreement CR 822038

### Field Operations Protocol

### EPA-Compendium

**Title:** Compendium of Methods for Analysis of Metals and VOCs in  
Water

**Source:** U.S. Environmental Protection Agency

U.S. Environmental Protection Agency  
Office of Research and Development  
Human Exposure & Atmospheric Sciences Division  
Human Exposure Research Branch

**Notice:** The U.S. Environmental Protection Agency (EPA), through its Office of Research and Development (ORD), partially funded and collaborated in the research described here. This protocol is part of the Quality Systems Implementation Plan (QSIP) that was reviewed by the EPA and approved for use in this demonstration/scoping study. Mention of trade names or commercial products does not constitute endorsement or recommendation by EPA for use.

COMPENDIUM OF METHODS FOR ANALYSIS OF METALS AND VOCs IN WATER

TABLE OF CONTENTS

Method 200.8	--	Determination of Trace Elements in Waters and Wastes by Inductively Coupled Plasma-Mass Spectrometry
Method 524.2	--	Measurement of Purgeable Organic Compounds in Water by Capillary Column Gas Chromatography/Mass Spectrometry
Arsenic Method	--	To be provided by EPA

CME  
8/25/92

**METHOD 524.2. MEASUREMENT OF PURGEABLE ORGANIC COMPOUNDS IN WATER BY  
CAPILLARY COLUMN GAS CHROMATOGRAPHY/MASS SPECTROMETRY**

**Revision 4.0**

**August 1992**

**A. Alford-Stevens, J. W. Eichelberger, and W. L. Budde  
Method 524, Revision 1.0 (1983)**

**R. W. Slater, Jr.  
Revision 2.0 (1986)**

**J. W. Eichelberger, and W. L. Budde  
Revision 3.0 (1989)**

**J. W. Eichelberger, J. W. Munch, and T. A. Bellar  
Revision 4.0 (1992)**

**ENVIRONMENTAL MONITORING SYSTEMS LABORATORY  
OFFICE OF RESEARCH AND DEVELOPMENT  
U.S. ENVIRONMENTAL PROTECTION AGENCY  
CINCINNATI, OHIO 45268**

## METHOD 524.2

### MEASUREMENT OF PURGEABLE ORGANIC COMPOUNDS IN WATER BY CAPILLARY COLUMN GAS CHROMATOGRAPHY/MASS SPECTROMETRY

#### 1. SCOPE AND APPLICATION

- 1.1 This is a general purpose method for the identification and simultaneous measurement of purgeable volatile organic compounds in surface water, ground water, and drinking water in any stage of treatment (1,2). The method is applicable to a wide range of organic compounds, including the four trihalomethane disinfection by-products, that have sufficiently high volatility and low water solubility to be removed from water samples with purge and trap procedures. The following compounds can be determined by this method.

<u>Compound</u>	<u>Chemical Abstract Service Registry Number</u>
Acetone*	67-64-1
Acrylonitrile*	107-13-1
Allyl chloride*	107-05-1
Benzene	71-43-2
Bromobenzene	108-86-1
Bromochloromethane	74-97-5
Bromodichloromethane	75-27-4
Bromoform	75-25-2
Bromomethane	74-83-9
2-Butanone*	78-93-3
n-Butylbenzene	104-51-8
sec-Butylbenzene	135-98-8
tert-Butylbenzene	98-06-6
Carbon disulfide*	75-15-0
Carbon tetrachloride	56-23-5
Chloroacetonitrile*	107-14-2
Chlorobenzene	108-90-7
1-Chlorobutane*	109-69-3
Chloroethane	75-00-3
Chloroform	67-66-3
Chloromethane	74-87-3
2-Chlorotoluene	95-49-8
4-Chlorotoluene	106-43-4
Dibromochloromethane	124-48-1
1,2-Dibromo-3-chloropropane	96-12-8
1,2-Dibromoethane	106-93-4
Dibromomethane	74-95-3
1,2-Dichlorobenzene	95-50-1
1,3-Dichlorobenzene	541-73-1
1,4-Dichlorobenzene	106-46-7
trans-1,4-Dichloro-2-butene*	110-57-6
Dichlorodifluoromethane	75-71-8

1,1-Dichloroethane	75-34-3
1,2-Dichloroethane	107-06-2
1,1-Dichloroethene	75-35-4
cis-1,2-Dichloroethene	156-59-4
trans-1,2-Dichloroethene	156-60-5
1,2-Dichloropropane	78-87-5
1,3-Dichloropropane	142-28-9
2,2-Dichloropropane	590-20-7
1,1-Dichloropropene	563-58-6
1,1-Dichloropropanone*	513-88-2
cis-1,3-Dichloropropene	10061-01-5
trans-1,3-Dichloropropene	10061-02-6
Diethyl ether*	60-29-7
Ethylbenzene	100-41-4
Ethyl methacrylate*	97-63-2
Hexachlorobutadiene	87-68-3
Hexachloroethane*	67-72-1
2-Hexanone*	591-78-6
Isopropylbenzene	98-82-8
4-Isopropyltoluene	99-87-6
Methacrylonitrile*	126-98-7
Methylacrylate*	96-33-3
Methylene chloride	75-09-2
Methyl iodide*	74-88-4
Methylmethacrylate*	80-62-6
4-Methyl-2-pentanone*	108-10-1
Methyl-t-butyl ether*	1634-04-4
Naphthalene	91-20-3
Nitrobenzene*	98-95-3
2-Nitropropane*	79-46-9
Pentachloroethane*	76-01-7
Propionitrile*	107-12-0
n-Propylbenzene	103-65-1
Styrene	100-42-5
1,1,1,2-Tetrachloroethane	100-42-5
1,1,2,2-Tetrachloroethane	630-20-6
Tetrachloroethene	79-34-5
Tetrahydrofuran*	127-18-4
Toluene	109-99-9
1,2,3-Trichlorobenzene	108-88-3
1,2,4-Trichlorobenzene	87-61-6
1,1,1-Trichloroethane	120-82-1
1,1,2-Trichloroethane	71-55-6
Trichloroethene	79-00-5
Trichlorofluoromethane	79-01-6
1,2,3-Trichloropropane	75-69-4
1,2,4-Trimethylbenzene	96-18-4
1,3,5-Trimethylbenzene	95-63-6
Vinyl chloride	108-67-8
o-Xylene	75-01-4
m-Xylene	95-47-6
p-Xylene	108-38-3
	106-42-3

\* New Compound in Revision 4.0

1.2 Method detection limits (MDLs) (3) are compound, instrument and especially matrix dependent and vary from approximately 0.02 to 1.6  $\mu\text{g/L}$ . The applicable concentration range of this method is primarily column and matrix dependent, and is approximately 0.02 to 200  $\mu\text{g/L}$  when a wide-bore thick-film capillary column is used. Narrow-bore thin-film columns may have a capacity which limits the range to about 0.02 to 20  $\mu\text{g/L}$ . Volatile water soluble, polar compounds which have relatively low purging efficiencies can be determined using this method. Such compounds may be more susceptible to matrix effects, and the quality of the data may be adversely influenced.

1.3 Analytes that are not separated chromatographically, but which have different mass spectra and noninterfering quantitation ions (Table 1), can be identified and measured in the same calibration mixture or water sample as long as their concentrations are somewhat similar (Sect. 11.6.2). Analytes that have very similar mass spectra cannot be individually identified and measured in the same calibration mixture or water sample unless they have different retention times (Sect. 11.6.3). Coeluting compounds with very similar mass spectra, typically many structural isomers, must be reported as an isomeric group or pair. Two of the three isomeric xylenes and two of the three dichlorobenzenes are examples of structural isomers that may not be resolved on the capillary column, and if not, must be reported as isomeric pairs. The more water soluble compounds (> 2% solubility) and compounds with boiling points above 200°C are purged from the water matrix with lower efficiencies. These analytes may be more susceptible to matrix effects.

## 2. SUMMARY OF METHOD

2.1 Volatile organic compounds and surrogates with low water solubility are extracted (purged) from the sample matrix by bubbling an inert gas through the aqueous sample. Purged sample components are trapped in a tube containing suitable sorbent materials. When purging is complete, the sorbent tube is heated and backflushed with helium to desorb the trapped sample components into a capillary gas chromatography (GC) column interfaced to a mass spectrometer (MS). The column is temperature programmed to facilitate the separation of the method analytes which are then detected with the MS. Compounds eluting from the GC column are identified by comparing their measured mass spectra and retention times to reference spectra and retention times in a data base. Reference spectra and retention times for analytes are obtained by the measurement of calibration standards under the same conditions used for samples. The concentration of each identified component is measured by relating the MS response of the quantitation ion produced by that compound to the MS response of the quantitation ion produced by a compound that is used as an internal standard. Surrogate analytes, whose concentrations are known in every sample, are measured with the same internal standard calibration procedure.

## 3. DEFINITIONS

3.1 INTERNAL STANDARD (IS) -- A pure analyte(s) added to a sample, extract, or standard solution in known amount(s) and used to measure

the relative responses of other method analytes and surrogates that are components of the same sample or solution. The internal standard must be an analyte that is not a sample component.

- 3.2 SURROGATE ANALYTE (SA) -- A pure analyte(s), which is extremely unlikely to be found in any sample, and which is added to a sample aliquot in known amount(s) before extraction or other processing and is measured with the same procedures used to measure other sample components. The purpose of the SA is to monitor method performance with each sample.
- 3.3 LABORATORY DUPLICATES (LD1 and LD2) -- Two aliquots of the same sample taken in the laboratory and analyzed separately with identical procedures. Analyses of LD1 and LD2 indicates precision associated with laboratory procedures, but not with sample collection, preservation, or storage procedures.
- 3.4 FIELD DUPLICATES (FD1 and FD2) -- Two separate samples collected at the same time and place under identical circumstances and treated exactly the same throughout field and laboratory procedures. Analyses of FD1 and FD2 give a measure of the precision associated with sample collection, preservation and storage, as well as with laboratory procedures.
- 3.5 LABORATORY REAGENT BLANK (LRB) -- An aliquot of reagent water or other blank matrix that is treated exactly as a sample including exposure to all glassware, equipment, solvents, reagents, internal standards, and surrogates that are used with other samples. The LRB is used to determine if method analytes or other interferences are present in the laboratory environment, the reagents, or the apparatus.
- 3.6 FIELD REAGENT BLANK (FRB) -- An aliquot of reagent water or other blank matrix that is placed in a sample container in the laboratory and treated as a sample in all respects, including shipment to the sampling site, exposure to sampling site conditions, storage, preservation, and all analytical procedures. The purpose of the FRB is to determine if method analytes or other interferences are present in the field environment.
- 3.7 LABORATORY PERFORMANCE CHECK SOLUTION (LPC) -- A solution of one or more compounds (analytes, surrogates, internal standard, or other test compounds) used to evaluate the performance of the instrument system with respect to a defined set of method criteria.
- 3.8 LABORATORY FORTIFIED BLANK (LFB) -- An aliquot of reagent water or other blank matrix to which known quantities of the method analytes are added in the laboratory. The LFB is analyzed exactly like a sample, and its purpose is to determine whether the methodology is in control, and whether the laboratory is capable of making accurate and precise measurements.
- 3.9 LABORATORY FORTIFIED SAMPLE MATRIX (LFM) -- An aliquot of an environmental sample to which known quantities of the method analytes are

added in the laboratory. The LFM is analyzed exactly like a sample, and its purpose is to determine whether the sample matrix contributes bias to the analytical results. The background concentrations of the analytes in the sample matrix must be determined in a separate aliquot and the measured values in the LFM corrected for background concentrations.

- 3.10 STOCK STANDARD SOLUTION (SSS) -- A concentrated solution containing one or more method analytes prepared in the laboratory using assayed reference materials or purchased from a reputable commercial source.
- 3.11 PRIMARY DILUTION STANDARD SOLUTION (PDS) -- A solution of several analytes prepared in the laboratory from stock standard solutions and diluted as needed to prepare calibration solutions and other needed analyte solutions.
- 3.12 CALIBRATION STANDARD (CAL) -- A solution prepared from the primary dilution standard solution or stock standard solutions and the internal standards and surrogate analytes. The CAL solutions are used to calibrate the instrument response with respect to analyte concentration.
- 3.13 QUALITY CONTROL SAMPLE (QCS) -- A solution of method analytes of known concentrations which is used to fortify an aliquot of LRB or sample matrix. The QCS is obtained from a source external to the laboratory and different from the source of calibration standards. It is used to check laboratory performance with externally prepared test materials.

#### 4. INTERFERENCES

- 4.1 During analysis, major contaminant sources are volatile materials in the laboratory and impurities in the inert purging gas and in the sorbent trap. The use of Teflon tubing, Teflon thread sealants, or flow controllers with rubber components in the purging device should be avoided since such materials out-gas organic compounds which will be concentrated in the trap during the purge operation. Analyses of laboratory reagent blanks provide information about the presence of contaminants. When potential interfering peaks are noted in laboratory reagent blanks, the analyst should change the purge gas source and regenerate the molecular sieve purge gas filter. Subtracting blank values from sample results is not permitted.
- 4.2 Interfering contamination may occur when a sample containing low concentrations of volatile organic compounds is analyzed immediately after a sample containing relatively high concentrations of volatile organic compounds. A preventive technique is between-sample rinsing of the purging apparatus and sample syringes with two portions of reagent water. After analysis of a sample containing high concentrations of volatile organic compounds, one or more laboratory reagent blanks should be analyzed to check for cross-contamination.
- 4.3 Special precautions must be taken to determine methylene chloride. The analytical and sample storage area should be isolated from all



atmospheric sources of methylene chloride, otherwise random background levels will result. Since methylene chloride will permeate Teflon tubing, all GC carrier gas lines and purge gas plumbing should be constructed of stainless steel or copper tubing. Laboratory worker's clothing should be cleaned frequently since clothing previously exposed to methylene chloride fumes during common liquid/liquid extraction procedures can contribute to sample contamination.

- 4.4 Traces of ketones, methylene chloride, and some other organic solvents can be present even in the highest purity methanol. This is another potential source of contamination, and should be assessed before standards are prepared in the methanol.

## 5. SAFETY

- 5.1 The toxicity or carcinogenicity of chemicals used in this method has not been precisely defined; each chemical should be treated as a potential health hazard, and exposure to these chemicals should be minimized. Each laboratory is responsible for maintaining awareness of OSHA regulations regarding safe handling of chemicals used in this method. Additional references to laboratory safety are available (4-6) for the information of the analyst.
- 5.2 The following method analytes have been tentatively classified as known or suspected human or mammalian carcinogens: benzene, carbon tetrachloride, 1,4-dichlorobenzene, 1,2-dichloroethane, hexachlorobutadiene, 1,1,2,2-tetrachloroethane, 1,1,2-trichloroethane, chloroform, 1,2-dibromoethane, tetrachloroethene, trichloroethene, and vinyl chloride. Pure standard materials and stock standard solutions of these compounds should be handled in a hood. A NIOSH/MESA approved toxic gas respirator should be worn when the analyst handles high concentrations of these toxic compounds.

## 6. EQUIPMENT AND SUPPLIES

- 6.1 **SAMPLE CONTAINERS** -- 40-mL to 120-mL screw cap vials each equipped with a Teflon faced silicone septum. Prior to use, wash vials and septa with detergent and rinse with tap and distilled water. Allow the vials and septa to air dry at room temperature, place in a 105°C oven for 1 hr, then remove and allow to cool in an area known to be free of organics.
- 6.2 **PURGE AND TRAP SYSTEM** -- The purge and trap system consists of three separate pieces of equipment: purging device, trap, and desorber. Systems are commercially available from several sources that meet all of the following specifications.
- 6.2.1 The all glass purging device (Figure 1) should be designed to accept 25-mL samples with a water column at least 5 cm deep. A smaller (5-mL) purging device is recommended if the GC/MS system has adequate sensitivity to obtain the method detection limits required. Gaseous volumes above the sample must be kept to a minimum (< 15 mL) to eliminate dead volume

effects. A glass frit should be installed at the base of the sample chamber so the purge gas passes through the water column as finely divided bubbles with a diameter of  $< 3 \text{ mm}$  at the origin. Needle spargers may be used, however, the purge gas must be introduced at a point about  $5 \text{ mm}$  from the base of the water column. The use of a moisture control device is recommended to prohibit much of the trapped water vapor from entering the GC/MS and eventually causing instrumental problems.

6.2.2 The trap (Figure 2) must be at least  $25 \text{ cm}$  long and have an inside diameter of at least  $0.105 \text{ in.}$  Starting from the inlet, the trap should contain  $1.0 \text{ cm}$  of methyl silicone coated packing and the following amounts of adsorbents:  $1/3$  of 2,6-diphenylene oxide polymer,  $1/3$  of silica gel, and  $1/3$  of coconut charcoal. If it is not necessary to determine dichlorodifluoromethane, the charcoal can be eliminated and the polymer increased to fill  $2/3$  of the trap. Before initial use, the trap should be conditioned overnight at  $180^\circ\text{C}$  by backflushing with an inert gas flow of at least  $20 \text{ mL/min.}$  Vent the trap effluent to the room, not to the analytical column. Prior to daily use, the trap should be conditioned for  $10 \text{ min}$  at  $180^\circ\text{C}$  with backflushing. The trap may be vented to the analytical column during daily conditioning; however, the column must be run through the temperature program prior to analysis of samples. The use of alternative sorbents is acceptable, depending on the particular set of target analytes or other problems encountered, but the new trap packing must meet all quality control criteria described in Sect. 9.

6.2.3 The use of the methyl silicone coated packing is recommended, but not mandatory. The packing serves a dual purpose of protecting the Tenax adsorbant from aerosols, and also of insuring that the Tenax is fully enclosed within the heated zone of the trap thus eliminating potential cold spots. Alternatively, silanized glass wool may be used as a spacer at the trap inlet.

6.2.4 The desorber (Figure 2) must be capable of rapidly heating the trap to  $180^\circ\text{C}$  either prior to or at the beginning of the flow of desorption gas. The polymer section of the trap should not be heated higher than  $200^\circ\text{C}$  or the life expectancy of the trap will decrease. Trap failure is characterized by a pressure drop in excess of  $3 \text{ lb/in}^2$  across the trap during purging or by poor bromoform sensitivities. The desorber design illustrated in Fig. 2 meets these criteria.

### 6.3 GAS CHROMATOGRAPHY/MASS SPECTROMETER/DATA SYSTEM (GC/MS/DS)

6.3.1 The GC must be capable of temperature programming and should be equipped with variable-constant differential flow controllers so that the column flow rate will remain constant throughout desorption and temperature program operation. If the column oven is to be cooled to  $10^\circ\text{C}$  or lower, a subambient oven controller will likely be required. If syringe

injections of 4-bromofluorobenzene (BFB) will be used, a split/splitless injection port is required.

- 6.3.2 Capillary GC Columns. Any gas chromatography column that meets the performance specifications of this method may be used (Sect. 10.2.4.1). Separations of the calibration mixture must be equivalent or better than those described in this method. Four useful columns have been evaluated, and observed compound retention times for these columns are listed in Table 2.

6.3.2.1 Column 1 -- 60 m x 0.75 mm ID VOCOL (Supelco, Inc.) glass wide-bore capillary with a 1.5  $\mu$ m film thickness.

Column 2 -- 30 m x 0.53 mm ID DB-624 (J&W Scientific, Inc.) fused silica capillary with a 3  $\mu$ m film thickness.

Column 3 -- 30 m x 0.32 mm ID DB-5 (J&W Scientific, Inc.) fused silica capillary with a 1  $\mu$ m film thickness.

Column 4 -- 75 m x 0.53 mm id DB-624 (J&W Scientific, Inc.) fused silica capillary with a 3  $\mu$ m film thickness.

- 6.3.3 Interfaces between the GC and MS. The interface used depends on the column selected and the gas flow rate.

6.3.3.1 The wide-bore columns 1, 2, and 4 have the capacity to accept the standard gas flows from the trap during thermal desorption, and chromatography can begin with the onset of thermal desorption. Depending on the pumping capacity of the MS, an additional interface between the end of the column and the MS may be required. An open split interface (7) or an all-glass jet separator is an acceptable interface. Any interface can be used if the performance specifications described in this method (Sect. 9 and 10) can be achieved. The end of the transfer line after the interface, or the end of the analytical column if no interface is used, should be placed within a few mm of the MS ion source.

6.3.3.2 When narrow bore column 3 is used, a cryogenic interface placed just in front of the column inlet is suggested. This interface condenses the desorbed sample components in a narrow band on an uncoated fused silica precolumn using liquid nitrogen cooling. When all analytes have been desorbed from the trap, the interface is rapidly heated to transfer them to the analytical column. The end of the analytical column should be placed within a few mm of the MS ion source. A potential problem with this interface is blockage of the interface by frozen

water from the trap. This condition will result in a major loss in sensitivity and chromatographic resolution.

6.3.4 The mass spectrometer must be capable of electron ionization at a nominal electron energy of 70 eV. The spectrometer must be capable of scanning from 35 to 260 amu with a complete scan cycle time (including scan overhead) of 2 sec or less. (Scan cycle time = Total MS data acquisition time in seconds divided by number of scans in the chromatogram.) The spectrometer must produce a mass spectrum that meets all criteria in Table 3 when 25 ng or less of 4-bromofluorobenzene (BFB) is introduced into the GC. An average spectrum across the BFB GC peak may be used to test instrument performance.

6.3.5 An interfaced data system is required to acquire, store, reduce, and output mass spectral data. The computer software should have the capability of processing stored GC/MS data by recognizing a GC peak within any given retention time window, comparing the mass spectra from the GC peak with spectral data in a user-created data base, and generating a list of tentatively identified compounds with their retention times and scan numbers. The software must allow integration of the ion abundance of any specific ion between specified time or scan number limits. The software should also allow calculation of response factors as defined in Sect. 10.2.6 (or construction of a linear or second order regression calibration curve), calculation of response factor statistics (mean and standard deviation), and calculation of concentrations of analytes using either the calibration curve or the equation in Sect. 12.

#### 6.4 SYRINGE AND SYRINGE VALVES

6.4.1 Two 5-mL or 25-mL glass hypodermic syringes with Luer-Lok tip (depending on sample volume used).

6.4.2 Three 2-way syringe valves with Luer ends.

6.4.3 Micro syringes - 10, 100  $\mu$ L.

6.4.4 Syringes - 0.5, 1.0, and 5-mL, gas tight with shut-off valve.

#### 6.5 MISCELLANEOUS

6.5.1 Standard solution storage containers -- 15-mL bottles with Teflon lined screw caps.

### 7. REAGENTS AND STANDARDS

#### 7.1 TRAP PACKING MATERIALS

7.1.1 2,6-Diphenylene oxide polymer, 60/80 mesh, chromatographic grade (Tenax GC or equivalent).

- 7.1.2 Methyl silicone packing (optional) -- OV-1 (3%) on Chromosorb W, 60/80 mesh, or equivalent.
- 7.1.3 Silica gel -- 35/60 mesh, Davison, grade 15 or equivalent.
- 7.1.4 Coconut charcoal -- Prepare from Barnebey Cheney, CA-580-26 lot #M-2649 by crushing through 26 mesh screen.

## 7.2 REAGENTS

- 7.2.1 Methanol -- Demonstrated to be free of analytes.
- 7.2.2 Reagent water -- Prepare reagent water by passing tap water through a filter bed containing about 0.5 kg of activated carbon, by using a water purification system, or by boiling distilled water for 15 min followed by a 1-h purge with inert gas while the water temperature is held at 90°C. Store in clean, narrow-mouth bottles with Teflon lined septa and screw caps.
- 7.2.3 Hydrochloric acid (1+1) -- Carefully add measured volume of conc. HCl to equal volume of reagent water.
- 7.2.4 Vinyl chloride -- Certified mixtures of vinyl chloride in nitrogen and pure vinyl chloride are available from several sources (for example, Matheson, Ideal Gas Products, and Scott Gases).
- 7.2.5 Ascorbic acid -- ACS reagent grade, granular.
- 7.2.6 Sodium thiosulfate -- ACS reagent grade, granular.

- 7.3 STOCK STANDARD SOLUTIONS -- These solutions may be purchased as certified solutions or prepared from pure standard materials using the following procedures. One of these solutions is required for every analyte of concern, every surrogate, and the internal standard. A useful working concentration is about 1-5 mg/mL.

- 7.3.1 Place about 9.8 mL of methanol into a 10-mL ground-glass stoppered volumetric flask. Allow the flask to stand, unstoppered, for about 10 min or until all alcohol-wetted surfaces have dried and weigh to the nearest 0.1 mg.
- 7.3.2 If the analyte is a liquid at room temperature, use a 100- $\mu$ L syringe and immediately add two or more drops of reference standard to the flask. Be sure that the reference standard falls directly into the alcohol without contacting the neck of the flask. If the analyte is a gas at room temperature, fill a 5-mL valved gas-tight syringe with the standard to the 5.0-mL mark, lower the needle to 5 mm above the methanol meniscus, and slowly inject the standard into the neck area of the flask. The gas will rapidly dissolve in the methanol.

- 7.3.3 Reweigh, dilute to volume, stopper, then mix by inverting the flask several times. Calculate the concentration in  $\mu\text{g}/\mu\text{L}$  from the net gain in weight. When compound purity is certified at 96% or greater, the weight can be used without correction to calculate the concentration of the stock standard.
- 7.3.4 Store stock standard solutions in 15-mL bottles equipped with Teflon lined screw caps. Methanol solutions of acrylonitrile, methyl iodide, and methyl acrylate are stable for only one week at 4°C. Methanol solutions prepared from other liquid analytes are stable for at least 4 weeks when stored at 4°C. Methanol solutions prepared from gaseous analytes are not stable for more than 1 week when stored at < 0°C; at room temperature, they must be discarded after 1 day.
- 7.4 PRIMARY DILUTION STANDARDS -- Use stock standard solutions to prepare primary dilution standard solutions that contain all the analytes of concern in methanol or other suitable solvent. The primary dilution standards should be prepared at concentrations that can be easily diluted to prepare aqueous calibration solutions that will bracket the working concentration range. Store the primary dilution standard solutions with minimal headspace and check frequently for signs of deterioration or evaporation, especially just before preparing calibration solutions. Storage times described for stock standard solutions in Sect. 7.3.4 also apply to primary dilution standard solutions.
- 7.5 FORTIFICATION SOLUTIONS FOR INTERNAL STANDARD AND SURROGATES
- 7.5.1 A solution containing the internal standard and the surrogate compounds is required to prepare laboratory reagent blanks (also used as a laboratory performance check solution), and to fortify each sample. Prepare a fortification solution containing fluorobenzene (internal standard), 1,2-dichlorobenzene- $d_4$  (surrogate), and BFB (surrogate) in methanol at concentrations of 5  $\mu\text{g}/\text{mL}$  of each (any appropriate concentration is acceptable). A 5- $\mu\text{L}$  aliquot of this solution added to a 25-mL water sample volume gives concentrations of 1  $\mu\text{g}/\text{L}$  of each. A 5- $\mu\text{L}$  aliquot of this solution added to a 5-mL water sample volume gives a concentration of 5  $\mu\text{g}/\text{L}$  of each. Additional internal standards and surrogate analytes are optional. Additional surrogate compounds should be similar in physical and chemical characteristics to the analytes of concern.
- 7.6 PREPARATION OF LABORATORY REAGENT BLANK (LRB) -- Fill a 25-mL (or 5-mL) syringe with reagent water and adjust to the mark (no air bubbles). Inject an appropriate volume of the fortification solution containing the internal standard and surrogates through the Luer Lok valve into the reagent water. Transfer the LRB to the purging device. See Sect. 11.1.2.

7.7 PREPARATION OF LABORATORY FORTIFIED BLANK -- Prepare this exactly like a calibration standard (Sect. 7.8). This is a calibration standard that is treated as a sample.

## 7.8 PREPARATION OF CALIBRATION STANDARDS

7.8.1 The number of calibration solutions (CALs) needed depends on the calibration range desired. A minimum of three CAL solutions is required to calibrate a range of a factor of 20 in concentration. For a factor of 50, use at least four standards, and for a factor of 100 at least five standards. One calibration standard should contain each analyte of concern at a concentration of 2-10 times the method detection limit (Tables 4, 5, and 7) for that compound. The other CAL standards should contain each analyte of concern at concentrations that define the range of the method. Every CAL solution contains the internal standard and the surrogate compounds at the same concentration (5 µg/L suggested for a 5-mL sample; 1 µg/L for a 25-mL sample).

7.8.2 To prepare a calibration standard, add an appropriate volume of a primary dilution standard containing all analytes of concern to an aliquot of acidified (pH 2) reagent water in a volumetric flask. Also add an appropriate volume of internal standard and surrogate compound solution from Sect. 7.5.1. Use a microsyringe and rapidly inject the methanol solutions into the expanded area of the filled volumetric flask. Remove the needle as quickly as possible after injection. Mix by inverting the flask three times only. Discard the contents contained in the neck of the flask. Aqueous standards are not stable in a volumetric flask and should be discarded after 1 hr unless transferred to a sample bottle and sealed immediately.

## 8. SAMPLE COLLECTION, PRESERVATION, AND STORAGE

### 8.1 SAMPLE COLLECTION, DECHLORINATION, AND PRESERVATION

8.1.1 Collect all samples in duplicate. If samples, such as finished drinking water or waste water, are suspected to contain residual chlorine, add about 25 mg of ascorbic acid per 40 mL of sample to the sample bottle before filling. If the residual chlorine is likely to be present > 5 mg/L, a determination of the amount of the chlorine may be necessary. Diethyl-p-phenylenediamine (DPD) test kits are commercially available to determine residual chlorine in the field. Add an additional 25 mg of ascorbic acid per each 5 mg/L of residual chlorine. If compounds boiling below 25°C are not to be determined, sodium thiosulfate may be used to reduce the residual chlorine. Fill sample bottles to overflowing, but take care not to flush out the rapidly dissolving ascorbic acid. No air bubbles should pass through the sample as the bottle is filled, or be trapped in the sample when the bottle is sealed. Adjust the pH of the duplicate samples to

< 2 by carefully adding two drops of 1:1 HCl for each 40 mL of sample. Seal the sample bottles, Teflon face down, and shake vigorously for 1 min. Do not mix the ascorbic acid or sodium thiosulfate with the HCl prior to sampling.

8.1.2 When sampling from a water tap, open the tap and allow the system to flush until the water temperature has stabilized (usually about 10 min). Adjust the flow to about 500 mL/min and collect duplicate samples from the flowing stream.

8.1.3 When sampling from an open body of water, such as surface water, waste water, and possible leachate samples, partially fill a 1-quart wide-mouth bottle or 1-L beaker with sample from a representative area. Fill a 60 mL or a 120 mL sample vial with sample from the larger container, and adjust the pH of the sample to about 2 by adding 1:1 HCl dropwise while stirring. Check the pH with narrow range (1.4 to 2.8) pH paper. Record the number of drops of acid necessary to adjust the pH to 2. To collect actual samples, refill the large container with fresh sample and pour sample into sample vials. Follow filling instructions in Sect. 8.1.1. Add the appropriate number of drops of 1:1 HCl to each sample to adjust the pH to about 2. If samples are suspected to contain residual chlorine, add ascorbic acid or sodium thiosulfate according to Sect. 8.1.1.

8.1.4 The samples must be chilled to about 4°C when collected and maintained at that temperature until analysis. Field samples that will not be received at the laboratory on the day of collection must be packaged for shipment with sufficient ice to ensure that they will arrive at the laboratory with a substantial amount of ice remaining in the cooler.

8.1.5 If a sample foams vigorously when HCl is added, discard that sample. Collect a set of duplicate samples but do not acidify them. These samples must be flagged as "not acidified" and must be stored at 4°C or below. These samples must be analyzed within 24 hr of collection time.

## 8.2 SAMPLE STORAGE

8.2.1 Store samples at  $\leq 4^{\circ}\text{C}$  until analysis. The sample storage area must be free of organic solvent vapors and direct or intense light.

8.2.2 Analyze all samples within 14 days of collection. Samples not analyzed within this period must be discarded and replaced.

## 8.3 FIELD REAGENT BLANKS (FRB)

8.3.1 Duplicate FRBs must be handled along with each sample set, which is composed of the samples collected from the same



general sample site at approximately the same time. At the laboratory, fill field blank sample bottles with reagent water and sample preservatives, seal, and ship to the sampling site along with empty sample bottles and back to the laboratory with filled sample bottles. Wherever a set of samples is shipped and stored, it is accompanied by appropriate blanks. FRBs must remain hermetically sealed until analysis.

- 8.3.2 Use the same procedures used for samples to add ascorbic acid and HCl to blanks (Sect. 8.1.1). The same batch of ascorbic acid and HCl should be used for the field reagent blanks in the field.

## **9. QUALITY CONTROL**

- 9.1 Quality control (QC) requirements are the initial demonstration of laboratory capability followed by regular analyses of laboratory reagent blanks, field reagent blanks, and laboratory fortified blanks. Each laboratory must maintain records to document the quality of the data generated. Additional quality control practices are recommended.
- 9.2 Initial demonstration of low system background. Before any samples are analyzed, it must be demonstrated that a laboratory reagent blank (LRB) is reasonably free of contamination that would prevent the determination of any analyte of concern. Sources of background contamination are glassware, purge gas, sorbents, and equipment. Background contamination must be reduced to an acceptable level before proceeding with the next section. In general, background from method analytes should be below the method detection limit.
- 9.3 Initial demonstration of laboratory accuracy and precision. Analyze five to seven replicates of a laboratory fortified blank containing each analyte of concern at a concentration in the range of 0.2-5 µg/L (see appropriate regulations and maximum contaminant levels for guidance on appropriate concentrations).
- 9.3.1 Prepare each replicate by adding an appropriate aliquot of a quality control sample to reagent water. If a quality control sample containing the method analytes is not available, a primary dilution standard made from a source of reagents different than those used to prepare the calibration standards may be used. Also add the appropriate amounts of internal standard and surrogate compounds. Analyze each replicate according to the procedures described in Sect. 11, and on a schedule that results in the analyses of all replicates over a period of several days.
- 9.3.2 Calculate the measured concentration of each analyte in each replicate, the mean concentration of each analyte in all replicates, and mean accuracy (as mean percentage of true

value) for each analyte, and the precision (as relative standard deviation, RSD) of the measurements for each analyte. Calculate the MDL of each analyte using the equation described in Sect. 13.2 (3).

9.3.3 For each analyte, the mean accuracy, expressed as a percentage of the true value, should be 80-120% and the RSD should be < 20%. Some analytes, particularly the early eluting gases and late eluting higher molecular weight compounds, are measured with less accuracy and precision than other analytes. The MDLs must be sufficient to detect analytes at the required levels according to the SDWA Regulations. If these criteria are not met for an analyte, take remedial action and repeat the measurements for that analyte to demonstrate acceptable performance before samples are analyzed.

9.3.4 Develop and maintain a system of control charts to plot the precision and accuracy of analyte and surrogate measurements as a function of time. Charting surrogate recoveries is an especially valuable activity because surrogates are present in every sample and the analytical results will form a significant record of data quality.

9.4 Monitor the integrated areas of the quantitation ions of the internal standards and surrogates (Table 1) in all samples, continuing calibration checks, and blanks. These should remain reasonably constant over time. An abrupt change may indicate a matrix effect or an instrument problem. If a cryogenic interface is utilized, it may indicate an inefficient transfer from the trap to the column. These samples must be reanalyzed or a laboratory fortified duplicate sample analyzed to test for matrix effect. A more gradual drift of more than 50% in any area is indicative of a loss in sensitivity, and the problem must be found and corrected.

9.5 LABORATORY REAGENT BLANKS (LRB) -- With each batch of samples processed as a group within a work shift, analyze a LRB to determine the background system contamination. A FRB (Sect. 9.7) may be used in place of a LRB.

9.6 With each batch of samples processed as a group within a work shift, analyze a single laboratory fortified blank (LFB) containing each analyte of concern at a concentration as determined in Sect. 9.3. If more than 20 samples are included in a batch, analyze one LFB for every 20 samples. Use the procedures described in Sect. 9.3.3 to evaluate the accuracy of the measurements, and to estimate whether the MDLs can be obtained. If acceptable accuracy and MDLs cannot be achieved, the problem must be located and corrected before further samples are analyzed. Add these results to the ongoing control charts to document data quality.

9.7 With each set of field samples a field reagent blank (FRB) should be analyzed. The results of these analyses will help define contamination.

tion resulting from field sampling and transportation activities. If the FRB shows unacceptable contamination, a LRB must be measured to define the source of the impurities.

- 9.8 At least quarterly, replicate LFBs should be analyzed to determine the precision of the laboratory measurements. Add these results to the ongoing control charts to document data quality.
- 9.9 At least quarterly, analyze a quality control sample (QCS) from an external source. If measured analyte concentrations are not of acceptable accuracy, check the entire analytical procedure to locate and correct the problem source.
- 9.10 Sample matrix effects have not been observed when this method is used with distilled water, reagent water, drinking water, or ground water. Therefore, analysis of a laboratory fortified sample matrix (LFM) is not required unless the criteria in Section 9.4 are not met. If matrix effects are observed or suspected to be causing low recoveries, analyze a laboratory fortified matrix sample for that matrix. The sample results should be flagged and the LFM results should be reported with them.
- 9.11 Numerous other quality control measures are incorporated into other parts of this procedure, and serve to alert the analyst to potential problems.

## **10. CALIBRATION AND STANDARDIZATION**

- 10.1 Demonstration and documentation of acceptable initial calibration is required before any samples are analyzed and is required intermittently throughout sample analysis as dictated by results of continuing calibration checks. After initial calibration is successful, a continuing calibration check is required at the beginning of each 8 hr. period during which analyses are performed. Additional periodic calibration checks are good laboratory practice.

### **10.2 INITIAL CALIBRATION**

- 10.2.1 Calibrate the mass and abundance scales of the MS with calibration compounds and procedures prescribed by the manufacturer with any modifications necessary to meet the requirements in Sect. 10.2.2.
- 10.2.2 Introduce into the GC (either by purging a laboratory reagent blank or making a syringe injection) 25 ng or less of BFB and acquire mass spectra for  $m/z$  35-260 at 70 eV (nominal). Use the purging procedure and/or GC conditions given in Sect. 11. If the spectrum does not meet all criteria in Table 3, the MS must be retuned and adjusted to meet all criteria before proceeding with calibration. An average spectrum across the GC peak may be used to evaluate the performance of the system.
- 10.2.3 Purge a medium CAL solution, (e.g., 10-20  $\mu\text{g/L}$ ) using the procedure given in Sect. 11.

10.2.4 Performance criteria for the medium calibration. Examine the stored GC/MS data with the data system software. Figures 3 and 4 shown acceptable total ion chromatograms.

10.2.4.1 GC performance. Good column performance will produce symmetrical peaks with minimum tailing for most compounds. If peaks are unusually broad, or if peaks are running together with little valleys between them, the wrong column has been selected or remedial action is probably necessary (Sect. 10.3.6).

10.2.4.2 MS sensitivity. The GC/MS/DS peak identification software should be able to recognize a GC peak in the appropriate retention time window for each of the compounds in calibration solution, and make correct tentative identifications. If fewer than 99% of the compounds are recognized, system maintenance is required. See Sect. 10.3.6.

10.2.5 If all performance criteria are met, purge an aliquot of each of the other CAL solutions using the same GC/MS conditions.

10.2.6 Calculate a response factor (RF) for each analyte and isomer pair for each CAL solution using the internal standard fluorobenzene. Table 1 contains suggested quantitation ions for all compounds. This calculation is supported in acceptable GC/MS data system software (Sect. 6.3.5), and many other software programs. RF is a unitless number, but units used to express quantities of analyte and internal standard must be equivalent.

$$RF = \frac{(A_x)(Q_{is})}{(A_{is})(Q_x)}$$

where:  $A_x$  = integrated abundance of the quantitation ion of the analyte.  
 $A_{is}$  = integrated abundance of the quantitation ion of the internal standard.  
 $Q_x$  = quantity of analyte purged in nanograms or concentration units.  
 $Q_{is}$  = quantity of internal standard purged in ng or concentration units.

10.2.6.1 For each analyte and surrogate, calculate the mean RF from analyses of CAL solutions. Calculate the standard deviation (SD) and the relative standard deviation (RSD) from each mean:  $RSD = 100 (SD/M)$ . If the RSD of any analyte or surrogate mean RF exceeds 20%, either analyze additional aliquots of appropriate CAL solutions to obtain an acceptable RSD of RFs over the entire concentration range, or

take action to improve GC/MS performance Sect. 10.3.6). Surrogate compounds are present at the same concentration on every sample, calibration standard, and all types of blanks.

10.2.7 As an alternative to calculating mean response factors and applying the RSD test, use the GC/MS data system software or other available software to generate a linear or second order regression calibration curve.

10.3 CONTINUING CALIBRATION CHECK -- Verify the MS tune and initial calibration at the beginning of each 8-hr work shift during which analyses are performed using the following procedure.

10.3.1 Introduce into the GC (either by purging a laboratory reagent blank or making a syringe injection) 25 ng or less of BFB and acquire a mass spectrum that includes data for  $m/z$  35-260. If the spectrum does not meet all criteria (Table 3), the MS must be retuned and adjusted to meet all criteria before proceeding with the continuing calibration check.

10.3.2 Purge a medium concentration CAL solution and analyze with the same conditions used during the initial calibration.

10.3.3 Demonstrate acceptable performance for the criteria shown in Sect. 10.2.4.

10.3.4 Determine that the absolute areas of the quantitation ions of the internal standard and surrogates have not decreased by more than 30% from the areas measured in the most recent continuing calibration check, or by more than 50% from the areas measured during initial calibration. If these areas have decreased by more than these amounts, adjustments must be made to restore system sensitivity. These adjustments may require cleaning of the MS ion source, or other maintenance as indicated in Sect. 10.3.6, and recalibration. Control charts are useful aids in documenting system sensitivity changes.

10.3.5 Calculate the RF for each analyte of concern and surrogate compound from the data measured in the continuing calibration check. The RF for each analyte and surrogate must be within 30% of the mean value measured in the initial calibration. Alternatively, if a linear or second order regression is used, the concentration measured using the calibration curve must be within 30% of the true value of the concentration in the medium calibration solution. If these conditions do not exist, remedial action must be taken which may require recalibration.

- 10.3.6 Some possible remedial actions. Major maintenance such as cleaning an ion source, cleaning quadrupole rods, etc. require returning to the initial calibration step.
- 10.3.6.1 Check and adjust GC and/or MS operating conditions; check the MS resolution, and calibrate the mass scale.
  - 10.3.6.2 Clean or replace the splitless injection liner; silanize a new injection liner. This applies only if the injection liner is an integral part of the system.
  - 10.3.6.3 Flush the GC column with solvent according to manufacturer's instructions.
  - 10.3.6.4 Break off a short portion (about 1 meter) of the column from the end near the injector; or replace GC column. This action will cause a slight change in retention times. Analyst may need to redefine retention windows.
  - 10.3.6.5 Prepare fresh CAL solutions, and repeat the initial calibration step.
  - 10.3.6.6 Clean the MS ion source and rods (if a quadrupole).
  - 10.3.6.7 Replace any components that allow analytes to come into contact with hot metal surfaces.
  - 10.3.6.8 Replace the MS electron multiplier, or any other faulty components.
  - 10.3.6.9 Replace the trap, especially when only a few compounds fail the criteria in Sect. 10.3.5 while the majority are determined successfully. Also check for gas leaks in the purge and trap unit as well as the rest of the analytical system.
- 10.4 Optional calibration for vinyl chloride using a certified gaseous mixture of vinyl chloride in nitrogen can be accomplished by the following steps.
- 10.4.1 Fill the purging device with 25.0 mL (or 5-mL) of reagent water or aqueous calibration standard.
  - 10.4.2 Start to purge the aqueous mixture. Inject a known volume (between 100 and 2000  $\mu$ L) of the calibration gas (at room temperature) directly into the purging device with a gas tight syringe. Slowly inject the gaseous sample through a septum seal at the top of the purging device at 2000  $\mu$ L/min. If the injection of the standard is made through the aqueous

sample inlet port, flush the dead volume with several mL of room air or carrier gas. Inject the gaseous standard before 5 min of the 11-min purge time have elapsed.

- 10.4.3 Determine the aqueous equivalent concentration of vinyl chloride standard, in  $\mu\text{g/L}$ , injected with the equation:

$$S = 0.102 (C)(V)$$

where  $S$  = Aqueous equivalent concentration of vinyl chloride standard in  $\mu\text{g/L}$ ;  
 $C$  = Concentration of gaseous standard in  $\text{mg/L (v/v)}$ ;  
 $V$  = Volume of standard injected in mL.

## 11. PROCEDURE

### 11.1 SAMPLE INTRODUCTION AND PURGING

- 11.1.1 This method is designed for a 25-mL sample volume, but a smaller (5 mL) sample volume is recommended if the GC/MS system has adequate sensitivity to achieve the required method detection limits. Adjust the helium purge gas flow rate to 40 mL/min. Attach the trap inlet to the purging device and open the syringe valve on the purging device.
- 11.1.2 Remove the plungers from two 25-mL (or 5-mL depending on sample size) syringes and attach a closed syringe valve to each. Warm the sample to room temperature, open the sample bottle, and carefully pour the sample into one of the syringe barrels to just short of overflowing. Replace the syringe plunger, invert the syringe, and compress the sample. Open the syringe valve and vent any residual air while adjusting the sample volume to 25.0-mL (or 5-mL). To all samples, blanks, and calibration standards, add 5- $\mu\text{L}$  (or an appropriate volume) of the fortification solution containing the internal standard and the surrogates to the sample through the syringe valve. Close the valve. Fill the second syringe in an identical manner from the same sample bottle. Reserve this second syringe for a reanalysis if necessary.
- 11.1.3 Attach the sample syringe valve to the syringe valve on the purging device. Be sure that the trap is cooler than 25°C, then open the sample syringe valve and inject the sample into the purging chamber. Close both valves and initiate purging. Purge the sample for 11.0 min at ambient temperature.
- 11.1.4 Standards and samples must be analyzed in exactly the same manner. Room temperature changes in excess of 10°F may adversely affect the accuracy and precision of the method.

## 11.2 SAMPLE DESORPTION

11.2.1 Non-cryogenic interface -- After the 11-min purge, place the purge and trap system in the desorb mode and preheat the trap to 180°C without a flow of desorption gas. Then simultaneously start the flow of desorption gas at a flow rate suitable for the column being used (optimum desorb flow rate is 15 mL/min) for about 4 min, begin the GC temperature program, and start data acquisition.

11.2.2 Cryogenic interface -- After the 11-min purge, place the purge and trap system in the desorb mode, make sure the cryogenic interface is a -150°C or lower, and rapidly heat the trap to 180°C while backflushing with an inert gas at 4 mL/min for about 5 min. At the end of the 5 min desorption cycle, rapidly heat the cryogenic trap to 250°C, and simultaneously begin the temperature program of the gas chromatograph, and start data acquisition.

11.2.3 While the trapped components are being introduced into the gas chromatograph (or cryogenic interface), empty the purging device using the sample syringe and wash the chamber with two 25-mL flushes of reagent water. After the purging device has been emptied, leave syringe valve open to allow the purge gas to vent through the sample introduction needle.

11.3 GAS CHROMATOGRAPHY/MASS SPECTROMETRY -- Acquire and store data over the nominal mass range 35-260 with a total cycle time (including scan overhead time) of 2 sec or less. If water, methanol, or carbon dioxide cause a background problem, start at 47 or 48 m/z. If ketones are to be determined, data must be acquired starting at m/z 43. Cycle time must be adjusted to measure five or more spectra during the elution of each GC peak. Suggested temperature programs are provided below. Alternative temperature programs can be used.

11.3.1 Single ramp linear temperature program for wide bore column 1 and 2 with a jet separator. Adjust the helium carrier gas flow rate to within the capacity of the separator, or about 15 mL/min. The column temperature is reduced 10°C and held for 5 min from the beginning of desorption, then programmed to 160°C at 6°C/min, and held until all components have eluted.

11.3.2 Multi-ramp temperature program for wide bore column 2 with the open split interface. Adjust the helium carrier gas flow rate to about 4.6 mL/min. The column temperature is reduced to 10°C and held for 6 min from the beginning of desorption, then heated to 70°C at 10°/min, heated to 120°C at 5°/min, heated to 180° at 8°/min, and held at 180° until all compounds have eluted.



- 11.3.3 Single ramp linear temperature program for narrow bore column 3 with a cryogenic interface. Adjust the helium carrier gas flow rate to about 4 mL/min. The column temperature is reduced to 10°C and held for 5 min from the beginning of vaporization from the cryogenic trap, programmed at 6°/min for 10 min, then 15°/min for 5 min to 145°C, and held until all components have eluted.
- 11.3.4 Multi-ramp temperature program for wide bore column 4 with the open split interface. Adjust the helium carrier gas flow rate to about 7.0 mL/min. The column temperature is - 10°C and held for 6 min. from beginning of desorption, then heated to 100°C at 10°C/min, heated to 200°C at 5°C/min and held at 200°C for 8 min or until all compounds of interest had eluted.
- 11.4 TRAP RECONDITIONING -- After desorbing the sample for 4 min, recondition the trap by returning the purge and trap system to the purge mode. Wait 15 sec, then close the syringe valve on the purging device to begin gas flow through the trap. Maintain the trap temperature at 180°C. Maintain the moisture control module, if utilized, at 90°C to remove residual water. After approximately 7 min, turn off the trap heater and open the syringe valve to stop the gas flow through the trap. When the trap is cool, the next sample can be analyzed.
- 11.5 TERMINATION OF DATA ACQUISITION -- When all the sample components have eluted from the GC, terminate MS data acquisition. Use appropriate data output software to display full range mass spectra and appropriate plots of ion abundance as a function of time. If any ion abundance exceeds the system working range, dilute the sample aliquot in the second syringe with reagent water and analyze the diluted aliquot.
- 11.6 IDENTIFICATION OF ANALYTES -- Identify a sample component by comparison of its mass spectrum (after background subtraction) to a reference spectrum in the user-created data base. The GC retention time of the sample component should be within three standard deviations of the mean retention time of the compound in the calibration mixture.
- 11.6.1 In general, all ions that are present above 10% relative abundance in the mass spectrum of the standard should be present in the mass spectrum of the sample component and should agree within absolute 20%. For example, if an ion has a relative abundance of 30% in the standard spectrum, its abundance in the sample spectrum should be in the range of 10 to 50%. Some ions, particularly the molecular ion, are of special importance, and should be evaluated even if they are below 10% relative abundance.
- 11.6.2 Identification requires expert judgment when sample components are not resolved chromatographically and produce mass spectra containing ions contributed by more than one analyte.

When GC peaks obviously represent more than one sample component (i.e., broadened peak with shoulder(s) or valley between two or more maxima), appropriate analyte spectra and background spectra can be selected by examining plots of characteristic ions for tentatively identified components. When analytes coelute (i.e., only one GC peak is apparent), the identification criteria can be met but each analyte spectrum will contain extraneous ions contributed by the coeluting compound. Because purgeable organic compounds are relatively small molecules and produce comparatively simple mass spectra, this is not a significant problem for most method analytes.

11.6.3 Structural isomers that produce very similar mass spectra can be explicitly identified only if they have sufficiently different GC retention times. Acceptable resolution is achieved if the height of the valley between two peaks is less than 25% of the average height of the two peaks. Otherwise, structural isomers are identified as isomeric pairs. Two of the three isomeric xylenes and two of the three dichlorobenzenes are examples of structural isomers that may not be resolved on the capillary columns. If unresolved, these groups of isomers must be reported as isomeric pairs.

11.6.4 Methylene chloride, acetone, carbon disulfide, and other background components appear in variable quantities in laboratory and field reagent blanks, and generally cannot be accurately measured. Subtraction of the concentration in the blank from the concentration in the sample is not acceptable because the concentration of the background in the blank is highly variable.

## 12. DATA ANALYSIS AND CALCULATIONS

12.1 Complete chromatographic resolution is not necessary for accurate and precise measurements of analyte concentrations if unique ions with adequate intensities are available for quantitation.

12.1.1 Calculate analyte and surrogate concentrations.

$$C_x = \frac{(A_x)(Q_{is}) 1000}{(A_{is}) RF V}$$

where:  $C_x$  = concentration of analyte or surrogate in  $\mu\text{g/L}$  in the water sample.  
 $A_x$  = integrated abundance of the quantitation ion of the analyte in the sample.  
 $A_{is}$  = integrated abundance of the quantitation ion of the internal standard in the sample.  
 $Q_{is}$  = total quantity (in micrograms) of internal standard added to the water sample.  
 $V$  = original water sample volume in mL.  
 $RF$  = mean response factor of analyte from the initial calibration.

12.1.2 Alternatively, use the GC/MS system software or other available proven software to compute the concentrations of the analytes and surrogates from the linear or second order regression curves.

12.1.3 Calculations should utilize all available digits of precision, but final reported concentrations should be rounded to an appropriate number of significant figures (one digit of uncertainty). Experience indicates that three significant figures may be used for concentrations above 99 µg/L, two significant figures for concentrations between 1- 99 µg/L, and one significant figure for lower concentrations.

12.1.4 Calculate the total trihalomethane concentration by summing the four individual trihalomethane concentrations.

### 13. METHOD PERFORMANCE

13.1 Single laboratory accuracy and precision data were obtained for the method analytes using laboratory fortified blanks with analytes at concentrations between 1 and 5 µg/L. Results were obtained using the four columns specified (Sect. 6.3.2.1) and the open split or jet separator (Sect. 6.3.3.1), or the cryogenic interface (Sect. 6.3.3.2). These data are shown in Tables 4-8.

13.2 With these data, method detection limits were calculated using the formula (3):

$$MDL = S t_{(n-1, 1-\alpha = 0.99)}$$

where:

$t_{(n-1, 1-\alpha = 0.99)}$  = Student's t value for the 99% confidence level with n-1 degrees of freedom,

n = number of replicates

S = the standard deviation of the replicate analyses.

### 14. POLLUTION PREVENTION

14.1 No solvents are utilized in this method except the extremely small volumes of methanol needed to make calibration standards. The only other chemicals used in this method are the neat materials in preparing standards and sample preservatives. All are used in extremely small amounts and pose no threat to the environment.

### 15. WASTE MANAGEMENT

15.1 There are no waste management issues involved with this method. Due to the nature of this method, the discarded samples are chemically less contaminated than when they were collected.

## 16. REFERENCES

1. A. Alford-Stevens, J.W. Eichelberger, W.L. Budde, "Purgeable Organic Compounds in Water by Gas Chromatography/Mass Spectrometry, Method 524." Environmental Monitoring and Support Laboratory, U.S. Environmental Protection Agency, Cincinnati, Ohio, February 1983.
2. C. Madding, "Volatile Organic Compounds in Water by Purge and Trap Capillary Column GC/MS," Proceedings of the Water Quality Technology Conference, American Water Works Association, Denver, CO, December 1984.
3. J.A. Glaser, D.L. Foerst, G.D. McKee, S.A. Quave, and W.L. Budde, "Trace Analyses for Wastewaters," Environ. Sci. Technol., **15**, 1426, 1981.
4. "Carcinogens-Working with Carcinogens," Department of Health, Education, and Welfare, Public Health Service, Center for Disease Control, National Institute for Occupational Safety and Health, Publication No. 77-206, August 1977.
5. "OSHA Safety and Health Standards, General Industry," (29CFR1910), Occupational Safety and Health Administration, OSHA 2206, (Revised, January 1976).
6. "Safety in Academic Chemistry Laboratories," American Chemical Society Publication, Committee on Chemical Safety, 3rd Edition, 1979.
7. R.F. Arrendale, R.F. Severson, and O.T. Chortyk, "Open Split Interface for Capillary Gas Chromatography/Mass Spectrometry," Anal. Chem. **56**, 1533, 1984.
8. J.J. Flesch, P.S. Fair, "The Analysis of Cyanogen Chloride in Drinking Water," Proceedings of Water Quality Technology Conference, American Water Works Association, St. Louis, MO., November 14-16, 1988.

# 17. TABLES, DIAGRAMS, FLOWCHARTS, AND VALIDATION DATA

TABLE 1. MOLECULAR WEIGHTS AND QUANTITATION IONS FOR METHOD ANALYTES

Compound	MW	Primary Quantitation Ion	Secondary Quantitation Ions
<u>Internal standard</u>			
Fluorobenzene	96	96	77
<u>Surrogates</u>			
4-Bromofluorobenzene	174	95	174,176
1,2-Dichlorobenzene-d4	150	152	115,150
<u>Target Analytes</u>			
Acetone	58	43	58
Acrylonitrile	53	52	53
Allyl chloride	76	76	49
Benzene	78	78	77
Bromobenzene	156	156	77,158
Bromochloromethane	128	128	49,130
Bromodichloromethane	162	83	85,127
Bromoform	250	173	175,252
Bromomethane	94	94	96
2-Butanone	72	43	57,72
n-Butylbenzene	134	91	134
sec-Butylbenzene	134	105	134
tert-Butylbenzene	134	119	91
Carbon disulfide	76	76	--
Carbon tetrachloride	152	117	119
Chloroacetonitrile	75	48	75
Chlorobenzene	112	112	77,114
1-Chlorobutane	92	56	49
Chloroethane	64	64	66
Chloroform	118	83	85
Chloromethane	50	50	52
2-Chlorotoluene	126	91	126
4-Chlorotoluene	126	91	126
Dibromochloromethane	206	129	127
1,2-Dibromo-3-Chloropropane	234	75	155,157
1,2-Dibromoethane	186	107	109,188
Dibromomethane	172	93	95,174
1,2-Dichlorobenzene	146	146	111,148
1,3-Dichlorobenzene	146	146	111,148
1,4-Dichlorobenzene	146	146	111,148

TABLE 1. (continued)

Compound	MW <sup>a</sup>	Primary Quantitation Ion	Secondary Quantitation Ions
trans-1,4-Dichloro-2-butene	124	53	88,75
Dichlorodifluoromethane	120	85	87
1,1-Dichloroethane	98	63	65,83
1,2-Dichloroethane	98	62	98
1,1-Dichloroethene	96	96	61,63
cis-1,2-Dichloroethene	96	96	61,98
trans-1,2-Dichloroethene	96	96	61,98
1,2-Dichloropropane	112	63	112
1,3-Dichloropropane	112	76	78
2,2-Dichloropropane	112	77	97
1,1-Dichloropropene	110	75	110,77
1,1-Dichloropropanone	126	43	83
cis-1,3-dichloropropene	110	75	110
trans-1,3-dichloropropene	110	75	110
Diethyl ether	74	59	45,73
Ethylbenzene	106	91	106
Ethyl methacrylate	114	69	99
Hexachlorobutadiene	258	225	260
Hexachloroethane	234	117	119,201
2-Hexanone	100	43	58
Isopropylbenzene	120	105	120
4-Isopropyltoluene	134	119	134,91
Methacrylonitrile	67	67	52
Methyl acrylate	86	55	85
Methylene chloride	84	84	86,49
Methyl iodide	142	142	127
Methylmethacrylate	100	69	99
4-Methyl-2-pentanone	100	43	58,85
Methyl-t-butyl ether	88	73	57
Naphthalene	128	128	--
Nitrobenzene	123	51	77
2-Nitropropane	89	46	--
Pentachloroethane	200	117	119,167
Propionitrile	55	54	--
n-Propylbenzene	120	91	120
Styrene	104	104	78
1,1,1,2-Tetrachloroethane	166	131	133,119
1,1,2,2-Tetrachloroethane	166	83	131,85
Tetrachloroethene	164	166	168,129
Tetrahydrofuran	72	71	72,42
Toluene	92	92	91
1,2,3-Trichlorobenzene	180	180	182
1,2,4-Trichlorobenzene	180	180	182
1,1,1-Trichloroethane	132	97	99,61
1,1,2-Trichloroethane	132	83	97,85

TABLE 1. (continued)

Compound	MW <sup>a</sup>	Primary Quantitation Ion	Secondary Quantitation Ions
Trichloroethene	130	95	130, 132
Trichlorofluoromethane	136	101	103
1,2,3-Trichloropropane	146	75	77
1,2,4-Trimethylbenzene	120	105	120
1,3,5-Trimethylbenzene	120	105	120
Vinyl Chloride	62	62	64
o-Xylene	106	106	91
m-Xylene	106	106	91
p-Xylene	106	106	91

<sup>a</sup>Monoisotopic molecular weight calculated from the atomic masses of the isotopes with the smallest masses.

TABLE 2. CHROMATOGRAPHIC RETENTION TIMES FOR METHOD ANALYTES  
ON THREE COLUMNS WITH FOUR SETS OF CONDITIONS<sup>a</sup>

Compound	Retention		Time (min:sec)		Column 4 <sup>e</sup>
	Column 1 <sup>b</sup>	Column 2 <sup>b</sup>	Column 2 <sup>c</sup>	Column 3 <sup>d</sup>	
<u>Internal standard</u>					
Fluorobenzene	8:49	6:27	14:06	8:03	22:00
<u>Surrogates</u>					
4-Bromofluorobenzene	18:38	15:43	23:38		31:21
1,2-Dichlorobenzene-d4	22:16	19:08	27:25		35:51
<u>Target Analytes</u>					
Acetone					16:14
Acrylonitrile					17:49
Allyl chloride					16:58
Benzene	8:14	5:40	13:30	7:25	21:32
Bromobenzene	18:57	15:52	24:00	16:25	31:52
Bromochloromethane	6:44	4:23	12:22	5:38	20:20
Bromodichloromethane	10:35	8:29	15:48	9:20	23:36
Bromoform	17:56	14:53	22:46	15:42	30:32
Bromomethane	2:01	0:58	4:48	1:17	12:26
2-Butanone					19:41
n-Butylbenzene	22:13	19:29	27:32	17:57	35:41
sec-Butylbenzene	20:47	18:05	26:08	17:28	34:04
tert-Butylbenzene	20:17	17:34	25:36	17:19	33:26
Carbon Disulfide					16:30
Carbon Tetrachloride	7:37	5:16	13:10	7:25	21:11
Chloroacetonitrile					23:51
Chlorobenzene	15:46	13:01	20:40	14:20	28:26
1-Chlorobutane					21:00
Chloroethane	2:05	1:01		1:27	
Chloroform	6:24	4:48	12:36	5:33	20:27
Chloromethane	1:38	0:44	3:24	0:58	9:11
2-Chlorotoluene	19:20	16:25	24:32	16:44	32:21
4-Chlorotoluene	19:30	16:43	24:46	16:49	32:38
Cyanogen chloride (8)				1:03	
Dibromochloromethane	14:23	11:51	19:12	12:48	26:57
1,2-Dibromo-3-Chloropropane	24:32	21:05		18:02	38:20
1,2-Dibromoethane	14:44	11:50	19:24	13:36	27:19
Dibromomethane	10:39	7:56	15:26	9:05	23:22
1,2-Dichlorobenzene	22:31	19:10	27:26	17:47	35:55
1,3-Dichlorobenzene	21:13	18:08	26:22	17:28	34:31
1,4-Dichlorobenzene	21:33	18:23	26:36	17:38	34:45
t-1,4-Dichloro-2-butene					31:44
Dichlorodifluoromethane	1:33	0:42	3:08	0:53	7:16
1,1-Dichloroethane	4:51	2:56	10:48	4:02	18:46



TABLE 2. (continued)

Compound	Retention		Time (min:sec)		
	Column 1 <sup>b</sup>	Column 2 <sup>b</sup>	Column 2 <sup>c</sup>	Column 3 <sup>d</sup>	Column 4 <sup>e</sup>
1,2-Dichloroethane	8:24	5:50	13:38	7:00	21:31
1,1-Dichloroethene	2:53	1:34	7:50	2:20	16:01
cis-1,2-Dichloroethene	6:11	3:54	11:56	5:04	19:53
trans-1,2-Dichloroethene	3:59	2:22	9:54	3:32	17:54
1,2-Dichloropropane	10:05	7:40	15:12	8:56	23:08
1,3-Dichloropropane	14:02	11:19	18:42	12:29	26:23
2,2-Dichloropropane	6:01	3:48	11:52	5:19	19:54
1,1-Dichloropropanone					24:52
1,1-Dichloropropene	7:49	5:17	13:06	7:10	21:08
cis-1,3-dichloropropene	11:58		16:42		24:24
trans-1,3-dichloropropene	13:46		17:54		25:33
Diethyl ether					15:31
Ethylbenzene	15:59	13:23	21:00	14:44	28:37
Ethyl Methacrylate					25:35
Hexachlorobutadiene	26:59	23:41	32:04	19:14	42:03
Hexachloroethane					36:45
Hexanone					26:23
Isopropylbenzene	18:04	15:28	23:18	16:25	30:52
4-Isopropyltoluene	21:12	18:31	26:30	17:38	34:27
Methacrylonitrile					20:15
Methylacrylate					20:02
Methylene Chloride	3:36	2:04	9:16	2:40	17:18
Methyl Iodide					16:21
Methylmethacrylate					23:08
4-Methyl-2-pentanone					24:38
Methyl-t-butyl ether					17:56
Naphthalene	27:10	23:31	32:12	19:04	42:29
Nitrobenzene					39:02
2-Nitropropane					23:58
Pentachloroethane					33:33
Propionitrile					19:58
n-Propylbenzene	19:04	16:25	24:20	16:49	32:00
Styrene	17:19	14:36	22:24	15:47	29:57
1,1,1,2-Tetrachloroethane	15:56	13:20	20:52	14:44	28:35
1,1,2,2-Tetrachloroethane	18:43	16:21	24:04	15:47	31:35
Tetrachloroethene	13:44	11:09	18:36	13:12	26:27
Tetrahydrofuran					20:26
Toluene	12:26	10:00	17:24	11:31	25:13
1,2,3-Trichlorobenzene	27:47	24:11	32:58	19:14	43:31
1,2,4-Trichlorobenzene	26:33	23:05	31:30	18:50	41:26
1,1,1-Trichloroethane	7:16	4:50	12:50	6:46	20:51
1,1,2-Trichloroethane	13:25	11:03	18:18	11:59	25:59
Trichloroethene	9:35	7:16	14:48	9:01	22:42
Trichlorofluoromethane	2:16	1:11	6:12	1:46	14:18
1,2,3-Trichloropropane	19:01	16:14	24:08	16:16	31:47
1,2,4-Trimethylbenzene	20:20	17:42	31:30	17:19	33:33

TABLE 2. (continued)

Compound	Retention		Time (min:sec)		
	Column 1 <sup>b</sup>	Column 2 <sup>b</sup>	Column 2 <sup>c</sup>	Column 3 <sup>d</sup>	Column 4 <sup>e</sup>
1,3,5-Trimethylbenzene	19:28	16:54	24:50	16:59	32:26
Vinyl chloride	1:43	0:47	3:56	1:02	10:22
o-Xylene	17:07	14:31	22:16	15:47	29:56
m-Xylene	16:10	13:41	21:22	15:18	28:53
p-Xylene	16:07	13:41	21:18	15:18	28:53

<sup>a</sup>Columns 1-4 are those given in Sect. 6.3.2.1; retention times were measured from the beginning of thermal desorption from the trap (columns 1-2, and 4) or from the beginning of thermal release from the cryogenic interface (column 3).

<sup>b</sup>GC conditions given in Sect. 11.3.1.

<sup>c</sup>GC conditions given in Sect. 11.3.2.

<sup>d</sup>GC conditions given in Sect. 11.3.3.

<sup>e</sup>GC conditions given in Sect. 11.3.4.

**TABLE 3. ION ABUNDANCE CRITERIA FOR 4-BROMOFLUOROBENZENE (BFB)**

<b>Mass (M/z)</b>	<b>Relative Abundance Criteria</b>
50	15 to 40% of mass 95
75	30 to 80% of mass 95
95	Base Peak, 100% Relative Abundance
96	5 to 9% of mass 95
173	< 2% of mass 174
174	> 50% of mass 95
175	5 to 9% of mass 174
176	> 95% but < 101% of mass 174
177	5 to 9% of mass 176

TABLE 4. ACCURACY AND PRECISION DATA FROM 16-31 DETERMINATIONS OF THE METHOD ANALYTES IN REAGENT WATER USING WIDE BORE CAPILLARY COLUMN 1<sup>a</sup>

Compound	True Conc. Range ( $\mu\text{g/L}$ )	Mean Accuracy (% of True Value)	Rel. Std. Dev. (%)	Method Det. Limit <sup>b</sup> ( $\mu\text{g/L}$ )
Benzene	0.1-10	97	5.7	0.04
Bromobenzene	0.1-10	100	5.5	0.03
Bromochloromethane	0.5-10	90	6.4	0.04
Bromodichloromethane	0.1-10	95	6.1	0.08
Bromoform	0.5-10	101	6.3	0.12
Bromomethane	0.5-10	95	8.2	0.11
n-Butylbenzene	0.5-10	100	7.6	0.11
sec-Butylbenzene	0.5-10	100	7.6	0.13
tert-Butylbenzene	0.5-10	102	7.3	0.14
Carbon tetrachloride	0.5-10	84	8.8	0.21
Chlorobenzene	0.1-10	98	5.9	0.04
Chloroethane	0.5-10	89	9.0	0.10
Chloroform	0.5-10	90	6.1	0.03
Chloromethane	0.5-10	93	8.9	0.13
2-Chlorotoluene	0.1-10	90	6.2	0.04
4-Chlorotoluene	0.1-10	99	8.3	0.06
Dibromochloromethane	0.1-10	92	7.0	0.05
1,2-Dibromo-3-chloropropane	0.5-10	83	19.9	0.26
1,2-Dibromoethane	0.5-10	102	3.9	0.06
Dibromomethane	0.5-10	100	5.6	0.24
1,2-Dichlorobenzene	0.1-10	93	6.2	0.03
1,3-Dichlorobenzene	0.5-10	99	6.9	0.12
1,4-Dichlorobenzene	0.2-20	103	6.4	0.03
Dichlorodifluoromethane	0.5-10	90	7.7	0.10
1,1-Dichloroethane	0.5-10	96	5.3	0.04
1,2-Dichloroethane	0.1-10	95	5.4	0.06
1,1-Dichloroethene	0.1-10	94	6.7	0.12
cis-1,2 Dichloroethene	0.5-10	101	6.7	0.12
trans-1,2-Dichloroethene	0.1-10	93	5.6	0.06
1,2-Dichloropropane	0.1-10	97	6.1	0.04
1,3-Dichloropropane	0.1-10	96	6.0	0.04
2,2-Dichloropropane	0.5-10	86	16.9	0.35
1,1-Dichloropropene	0.5-10	98	8.9	0.10
cis-1,2-Dichloropropene				
trans-1,2-Dichloropropene				
Ethylbenzene	0.1-10	99	8.6	0.06
Hexachlorobutadiene	0.5-10	100	6.8	0.11
Isopropylbenzene	0.5-10	101	7.6	0.15
4-Isopropyltoluene	0.1-10	99	6.7	0.12
Methylene chloride	0.1-10	95	5.3	0.03
Naphthalene	0.1-100	104	8.2	0.04
n-Propylbenzene	0.1-10	100	5.8	0.04
Styrene	0.1-100	102	7.2	0.04

TABLE 4. (Continued)

Compound	True Conc. Range ( $\mu\text{g/L}$ )	Mean Accuracy (% of True Value)	Rel. Std. Dev. (%)	Method Det. Limit <sup>b</sup> ( $\mu\text{g/L}$ )
1,1,1,2-Tetrachloroethane	0.5-10	90	6.8	0.05
1,1,2,2-Tetrachloroethane	0.1-10	91	6.3	0.04
Tetrachloroethene	0.5-10	89	6.8	0.14
Toluene	0.5-10	102	8.0	0.11
1,2,3-Trichlorobenzene	0.5-10	109	8.6	0.03
1,2,4-Trichlorobenzene	0.5-10	108	8.3	0.04
1,1,1-Trichloroethane	0.5-10	98	8.1	0.08
1,1,2-Trichloroethane	0.5-10	104	7.3	0.10
Trichloroethene	0.5-10	90	7.3	0.19
Trichlorofluoromethane	0.5-10	89	8.1	0.08
1,2,3-Trichloropropane	0.5-10	108	14.4	0.32
1,2,4-Trimethylbenzene	0.5-10	99	8.1	0.13
1,3,5-Trimethylbenzene	0.5-10	92	7.4	0.05
Vinyl chloride	0.5-10	98	6.7	0.17
o-Xylene	0.1-31	103	7.2	0.11
m-Xylene	0.1-10	97	6.5	0.05
p-Xylene	0.5-10	104	7.7	0.13

<sup>a</sup>Data obtained by using column 1 with a jet separator interface and a quadrupole mass spectrometer (Sect. 11.3.1) with analytes divided among three solutions.

<sup>b</sup>Replicate samples at the lowest concentration listed in column 2 of this table were analyzed. These results were used to calculate MDLs.

**TABLE 5. ACCURACY AND PRECISION DATA FROM SEVEN DETERMINATIONS OF METHOD ANALYTES IN REAGENT WATER USING THE CRYOGENIC TRAPPING OPTION AND A NARROW BORE CAPILLARY COLUMN 3°**

Compound	True Conc. (µg/L)	Mean Accuracy (% of True Value)	Rel. Std. Dev. (%)	Method Dect. Limit (µg/L)
Benzene	0.1	99	6.2	0.03
Bromobenzene	0.5	97	7.4	0.11
Bromochloromethane	0.5	97	5.8	0.07
Bromodichloromethane	0.1	100	4.6	0.03
Bromoform	0.1	99	5.4	0.20
Bromomethane	0.1	99	7.1	0.06
n-Butylbenzene	0.5	94	6.0	0.03
sec-Butylbenzene	0.5	90	7.1	0.12
tert-Butylbenzene	0.5	90	2.5	0.33
Carbon tetrachloride	0.1	92	6.8	0.08
Chlorobenzene	0.1	91	5.8	0.03
Chloroethane	0.1	100	5.8	0.02
Chloroform	0.1	95	3.2	0.02
Chloromethane	0.1	99	4.7	0.05
2-Chlorotoluene	0.1	99	4.6	0.05
4-Chlorotoluene	0.1	96	7.0	0.05
Cyanogen chloride <sup>b</sup>		92	10.6	0.30
Dibromochloromethane	0.1	99	5.6	0.07
1,2-Dibromo-3-chloropropane	0.1	92	10.0	0.05
1,2-Dibromoethane	0.1	97	5.6	0.02
Dibromomethane	0.1	93	6.9	0.03
1,2-Dichlorobenzene	0.1	97	3.5	0.05
1,3-Dichlorobenzene	0.1	99	6.0	0.05
1,4-Dichlorobenzene	0.1	93	5.7	0.04
Dichlorodifluoromethane	0.1	99	8.8	0.11
1,1-Dichloroethane	0.1	98	6.2	0.03
1,2-Dichloroethane	0.1	100	6.3	0.02
1,1-Dichloroethene	0.1	95	9.0	0.05
cis-1,2 Dichloroethene	0.1	100	3.7	0.06
trans-1,2-Dichloroethene	0.1	98	7.2	0.03
1,2-Dichloropropane	0.1	96	6.0	0.02
1,3-Dichloropropane	0.1	99	5.8	0.04
2,2-Dichloropropane	0.1	99	4.9	0.05
1,1-Dichloropropene	0.1	98	7.4	0.02
cis-1,3-Dichloropropene				
trans-1,3-Dichloropropene				
Ethylbenzene	0.1	99	5.2	0.03
Hexachlorobutadiene	0.1	100	6.7	0.04
Isopropylbenzene	0.5	98	6.4	0.10
4-Isopropyltoluene	0.5	87	13.0	0.26
Methylene chloride	0.5	97	13.0	0.09
Naphthalene	0.1	98	7.2	0.04

TABLE 5. (Continued)

Compound	True Conc. ( $\mu\text{g/L}$ )	Mean Accuracy (% of True Value)	Rel. Std. Dev. (%)	Method Dect. Limit ( $\mu\text{g/L}$ )
n-Propylbenzene	0.1	99	6.6	0.06
Styrene	0.1	96	19.0	0.06
1,1,1,2-Tetrachloroethane	0.1	100	4.7	0.04
1,1,2,2-Tetrachloroethane	0.5	100	12.0	0.20
Tetrachloroethene	0.1	96	5.0	0.05
Toluene	0.1	100	5.9	0.08
1,2,3-Trichlorobenzene	0.1	98	8.9	0.04
1,2,4-Trichlorobenzene	0.1	91	16.0	0.20
1,1,1-Trichloroethane	0.1	100	4.0	0.04
1,1,2-Trichloroethane	0.1	98	4.9	0.03
Trichloroethene	0.1	96	2.0	0.02
Trichlorofluoromethane	0.1	97	4.6	0.07
1,2,3-Trichloropropane	0.1	96	6.5	0.03
1,2,4-Trimethylbenzene	0.1	96	6.5	0.04
1,3,5-Trimethylbenzene	0.1	99	4.2	0.02
Vinyl chloride	0.1	96	0.2	0.04
o-Xylene	0.1	94	7.5	0.06
m-Xylene	0.1	94	4.6	0.03
p-Xylene	0.1	97	6.1	0.06

<sup>a</sup>Data obtained by using column 3 with a cryogenic interface and a quadrupole mass spectrometer (Sect 11.3.3).

<sup>b</sup>Reference 8.

TABLE 6. ACCURACY AND PRECISION DATA FROM SEVEN DETERMINATIONS  
OF THE METHOD ANALYTES IN REAGENT WATER USING WIDE BORE  
CAPILLARY COLUMN 2<sup>a</sup>

Compound	No. <sup>b</sup>	Mean Accuracy	RSD	Mean Accuracy	RSD
		(% of True Value, 2 µg/L Conc.)		(% of True Value, 0.2 µg/L Conc.)	
<u>Internal Standard</u>					
Fluorobenzene	1	-	-	-	-
<u>Surrogates</u>					
4-Bromofluorobenzene	2	98	1.8	96	1.3
1,2-Dichlorobenzene-d <sub>4</sub>	3	97	3.2	95	1.7
<u>Target Analytes</u>					
Benzene	37	97	4.4	113	1.8
Bromobenzene	38	102	3.0	101	1.9
Bromochloromethane	4	99	5.2	102	2.9
Bromodichloromethane	5	96	1.8	100	1.8
Bromoform	6	89	2.4	90	2.2
Bromomethane	7	55	27.	52	6.7
n-Butylbenzene	39	89	4.8	87	2.3
sec-Butylbenzene	40	102	3.5	100	2.8
tert-Butylbenzene	41	101	4.5	100	2.9
Carbon tetrachloride	8	84	3.2	92	2.6
Chlorobenzene	42	104	3.1	103	1.6
Chloroethane <sup>c</sup>					
Chloroform	9	97	2.0	95	2.1
Chloromethane	10	110	5.0	d	
2-Chlorotoluene	43	91	2.4	108	3.1
4-Chlorotoluene	44	89	2.0	108	4.4
Dibromochloromethane	11	95	2.7	100	3.0
1,2-Dibromo-3-chloropropane <sup>c</sup>					
1,2-Dibromoethane <sup>c</sup>					
Dibromomethane	13	99	2.1	95	2.2
1,2-Dichlorobenzene	45	93	2.7	94	5.1
1,3-Dichlorobenzene	46	100	4.0	87	2.3
1,4-Dichlorobenzene	47	98	4.1	94	2.8
Dichlorodifluoromethane	14	38	25.	d	
1,1-Dichloroethane	15	97	2.3	85	3.6
1,2-Dichloroethane	16	102	3.8	100	2.1
1,1-Dichloroethene	17	90	2.2	87	3.8
cis-1,2-Dichloroethene	18	100	3.4	89	2.9
trans-1,2-Dichloroethene	19	92	2.1	85	2.3



TABLE 6. (Continued)

Compound	No. <sup>b</sup>	Mean Accuracy (% of True Value, 2 µg/L Conc.)		RSD (%)	Mean Accuracy (% of True Value, 0.2 µg/L Conc.)		RSD (%)
1,2-Dichloropropane	20	102		2.2	103		2.9
1,3-Dichloropropane	21	92		3.7	93		3.2
2,2-Dichloropropane <sup>c</sup>							
1,1-Dichloropropene <sup>c</sup>							
cis-1,3-Dichloropropene <sup>c</sup>							
trans-1,3-Dichloropropene	25	96		1.7	99		2.1
Ethylbenzene	48	96		9.1	100		4.0
Hexachlorobutadiene	26	91		5.3	88		2.4
Isopropylbenzene	49	103		3.2	101		2.1
4-Isopropyltoluene	50	95		3.6	95		3.1
Methylene chloride	27	e			e		
Naphthalene	51	93		7.6	78		8.3
n-Propylbenzene	52	102		4.9	97		2.1
Styrene	53	95		4.4	104		3.1
1,1,1,2-Tetrachloroethane	28	99		2.7	95		3.8
1,1,2,2-Tetrachloroethane	29	101		4.6	84		3.6
Tetrachloroethene	30	97		4.5	92		3.3
Toluene	54	105		2.8	126		1.7
1,2,3-Trichlorobenzene	55	90		5.7	78		2.9
1,2,4-Trichlorobenzene	56	92		5.2	83		5.9
1,1,1-Trichloroethane	31	94		3.9	94		2.5
1,1,2-Trichloroethane	32	107		3.4	109		2.8
Trichloroethene	33	99		2.9	106		2.5
Trichlorofluoromethane	34	81		4.6	48		13.
1,2,3-Trichloropropane	35	97		3.9	91		2.8
1,2,4-Trimethylbenzene	57	93		3.1	106		2.2
1,3,5-Trimethylbenzene	58	88		2.4	97		3.2
Vinyl chloride	36	104		3.5	115		14.
o-Xylene	59	97		1.8	98		1.7
m-Xylene	60	f			f		
p-Xylene	61	98		2.3	103		1.4

<sup>a</sup>Data obtained using column 2 with the open split interface and an ion trap mass spectrometer (Sect. 11.3.2) with all method analytes in the same reagent water solution.

<sup>b</sup>Designation in Figures 1 and 2.

<sup>c</sup>Not measured; authentic standards were not available.

<sup>d</sup>Not found at 0.2 µg/L.

<sup>e</sup>Not measured; methylene chloride was in the laboratory reagent blank.

<sup>f</sup>m-xylene coelutes with and cannot be distinguished from its isomer p-xylene, No 61.

**TABLE 7. ACCURACY AND PRECISION DATA FROM SEVEN DETERMINATIONS  
OF METHOD ANALYTES IN REAGENT WATER USING WIDE BORE  
CAPILLARY COLUMN 4**

Compound	True Conc. ( $\mu\text{g/L}$ )	Mean Conc. Detected ( $\mu\text{g/L}$ )	Rel. Std. Dev. (%)	Method Det. Limit ( $\mu\text{g/L}$ )
Acetone	1.0	1.6	5.7%	0.28
Acrylonitrile	1.0	0.81	8.7%	0.22
Allyl chloride	1.0	0.90	4.7%	0.13
2-Butanone	2.0	2.7	5.6%	0.48
Carbon disulfide	0.20	0.19	15%	0.093
Chloroacetonitrile	1.0	0.83	4.7%	0.12
1-Chlorobutane	1.0	0.87	6.6%	0.18
t-1,2-Dichloro-2-butene	1.0	1.3	8.7%	0.36
1,1-Dichloropropanone	5.0	4.2	7.7%	1.0
Diethyl ether	1.0	0.92	9.5%	0.28
Ethyl methacrylate	0.20	0.23	3.9%	0.028
Hexachloroethane	0.20	0.18	10%	0.057
2-Hexanone	1.0	1.1	12%	0.39
Methacrylonitrile	1.0	0.92	4.2%	0.12
Methylacrylate	1.0	1.2	12%	0.45
Methyl iodide	0.20	0.19	3.1%	0.019
Methylmethacrylate	1.0	1.0	13%	0.43
4-Methyl-2-pentanone	0.40	0.56	9.7%	0.17
Methyl-tert-butylether	0.40	0.52	5.6%	0.090
Nitrobenzene	2.0	2.1	18%	1.2
2-Nitropropane	1.0	0.83	6.2%	0.16
Pentachloroethane	0.20	0.23	20%	0.14
Propionitrile	1.0	0.87	5.3%	0.14
Tetrahydrofuran	5.0	3.9	13%	1.6

**TABLE 8. ACCURACY AND PRECISION FROM FOUR DETERMINATIONS OF METHOD ANALYTES IN THREE WATER MATRICES FORTIFIED AT 20 µg/L**

Compound	REAGENT WATER			RAW WATER			TAP WATER		
	Mean (µg/L)	Dev. (%)	(% of True Value)	Mean (µg/L)	Dev. (%)	(% of True Value)	Mean (µg/L)	Dev. (%)	(% of True Value)
Acetone	19	12%	95%	21	3.7%	105%	22	8.2%	110%
Acrylonitrile	20	4.7%	100%	22	3.4%	110%	21	1.3%	105%
Allyl chloride	20	5.1%	100%	20	2.8%	100%	19	3.5%	95%
2-Butanone	17	11%	85%	19	7.3%	95%	17	5.6%	85%
Carbon disulfide	19	6.4%	95%	18	2.5%	90%	18	3.0%	90%
Chloroacetonitrile	20	4.1%	100%	23	4.7%	115%	23	1.3%	115%
1-Chlorobutane	18	6.4%	90%	19	2.2%	95%	17	2.2%	85%
t-1,2-Dichloro-2-butene	19	4.1%	95%	22	2.9%	110%	21	0.90%	105%
1,1-Dichloropropanone	20	5.6%	100%	22	6.4%	110%	21	7.7%	105%
Diethyl ether	18	6.7%	90%	22	3.4%	110%	22	2.6%	110%
Ethyl methacrylate	20	3.7%	100%	23	2.6%	115%	22	1.8%	110%
Hexachloroethane	20	6.1%	100%	21	2.5%	105%	21	2.0%	105%
2-Hexanone	19	6.3%	95%	21	3.8%	105%	21	4.0%	105%
Methacrylonitrile	20	3.4%	100%	23	2.9%	115%	22	2.0%	110%
Methylacrylate	20	3.7%	100%	22	3.1%	110%	21	2.1%	105%

TABLE 8 (Continued)

Compound	REAGENT WATER			RAW WATER			TAP WATER		
	Mean ( $\mu\text{g/L}$ )	Dev. (%)	(% of True Value)	Mean ( $\mu\text{g/L}$ )	Dev. (%)	(% of True Value)	Mean ( $\mu\text{g/L}$ )	Dev. (%)	(% of True Value)
Methyl iodide	20	4.4%	100%	19	3.8%	95%	19	3.0%	95%
Methylmethacrylate	20	3.7%	100%	23	3.3%	115%	23	2.7%	115%
4-Methyl-2-pentanone	19	8.7%	95%	21	5.5%	105%	22	7.2%	110%
Methyl-tert-butylether	19	3.5%	95%	22	2.5%	110%	22	3.6%	110%
Nitrobenzene	20	5.4%	100%	22	4.8%	110%	21	2.4%	105%
2-Nitropropane	20	6.1%	100%	23	5.1%	115%	22	3.2%	110%
Pentachloroethane	19	5.2%	95%	21	2.6%	105%	22	1.7%	110%
Propionitrile	20	4.5%	100%	23	3.9%	115%	23	2.4%	115%
Tetrahydrofuran	20	2.8%	100%	24	3.2%	120%	21	2.9%	105%

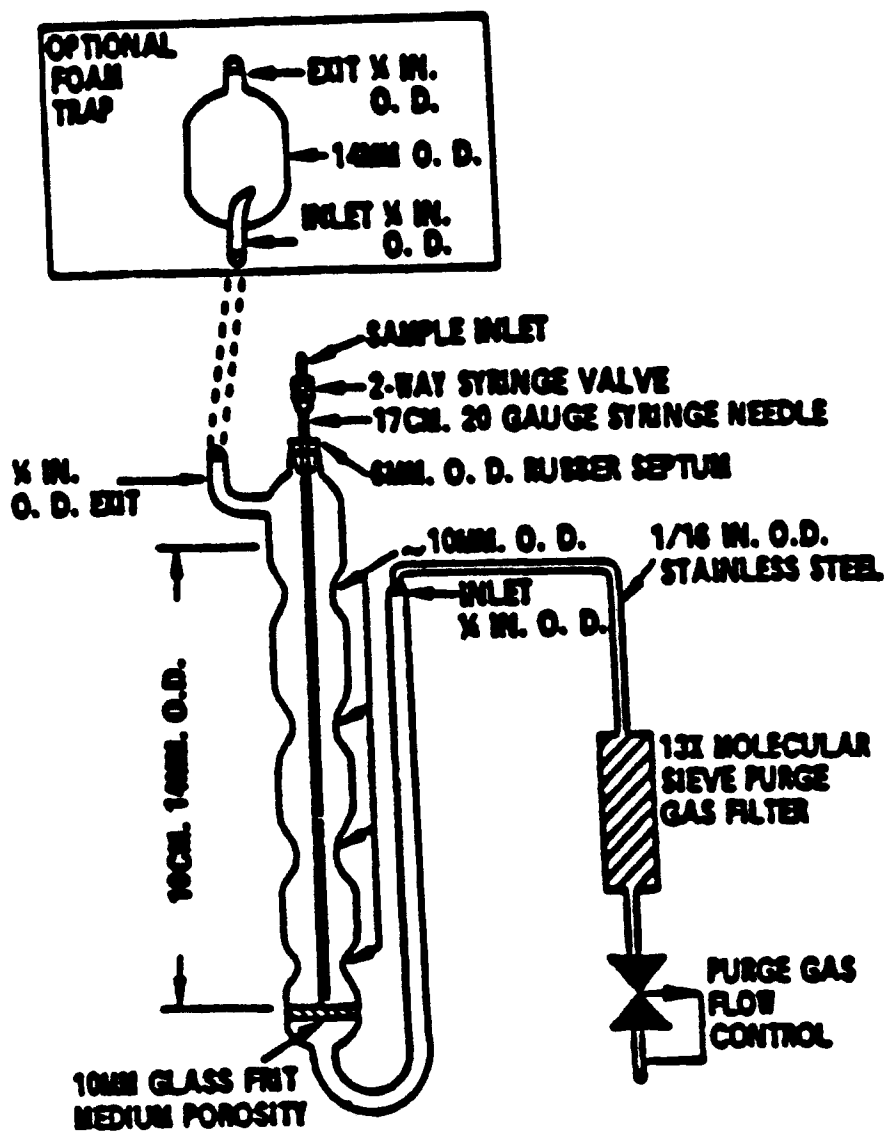


FIGURE 1. PURGING DEVICE



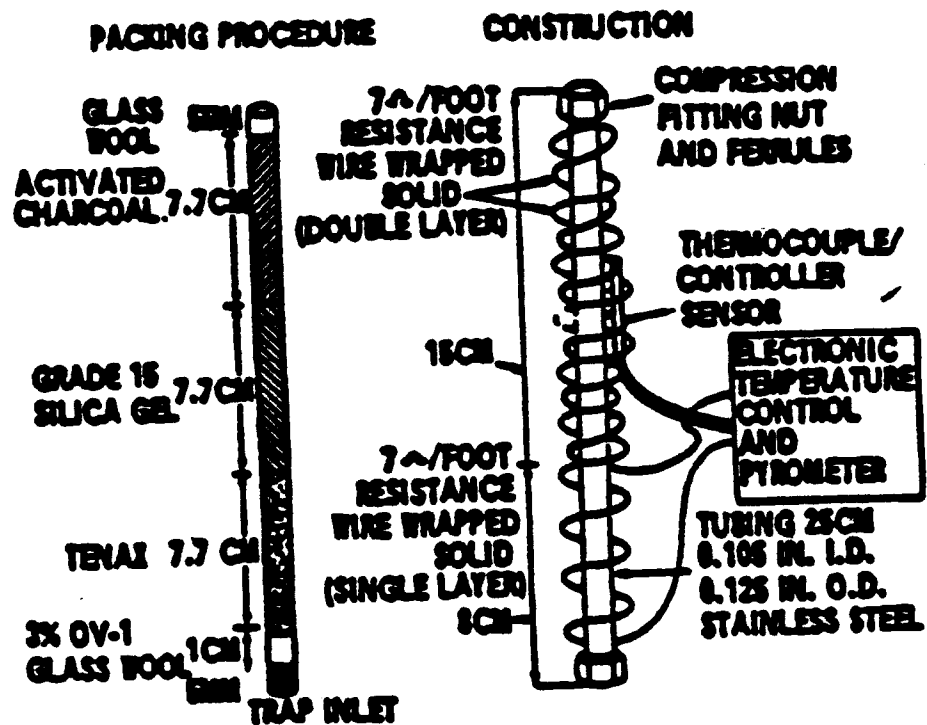


FIGURE 2. TRAP PACKINGS AND CONSTRUCTION TO INCLUDE DESORB CAPABILITY

FIGURE 3. NORMALIZED TOTAL ION CURRENT CHROMATOGRAM FROM A VOLATILE COMPOUND CALIBRATION MIXTURE CONTAINING 25  $\mu$ g (5  $\mu$ g/L) OF MOST COMPOUNDS. THE COMPOUND IDENTIFICATION NUMBERS ARE GIVEN IN TABLE 6.



**FIGURE 4. AMPLIFIED FIRST EIGHT MINUTES OF A TOTAL ION CURRENT CHROMATOGRAM FROM A VOLATILE COMPOUND CALIBRATION MIXTURE CONTAINING 25  $\mu$ g (5  $\mu$ g/L) OF EACH COMPONENT. THE COMPOUND IDENTIFICATION NUMBERS ARE GIVEN IN TABLE 6.**

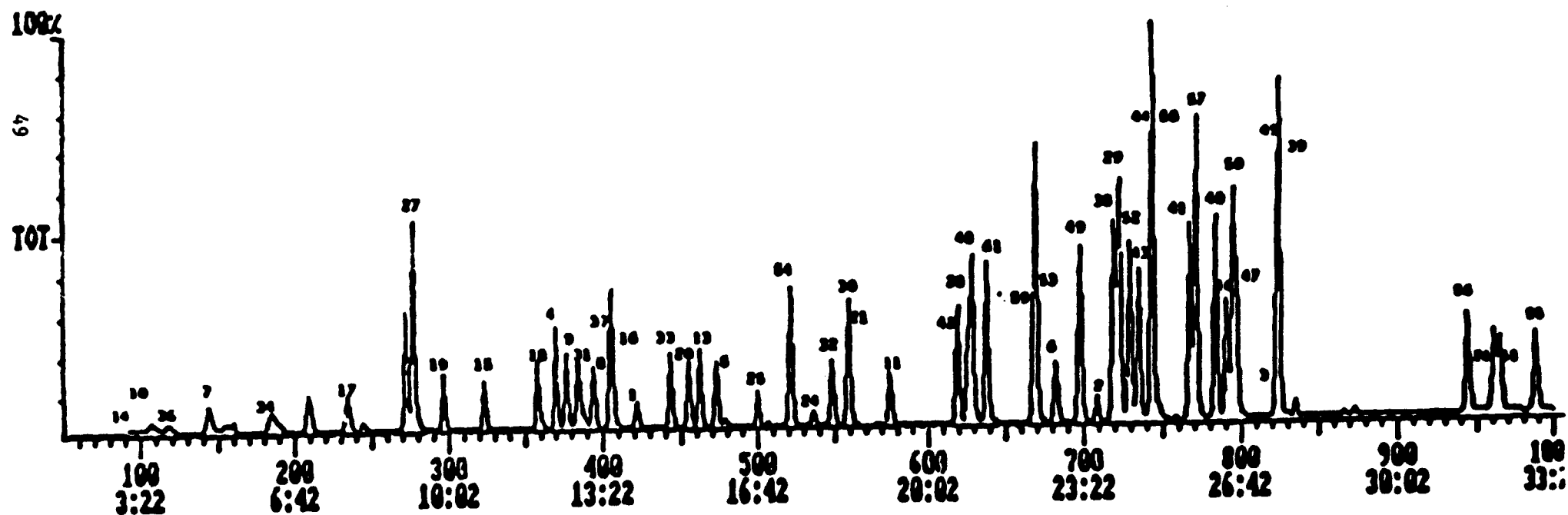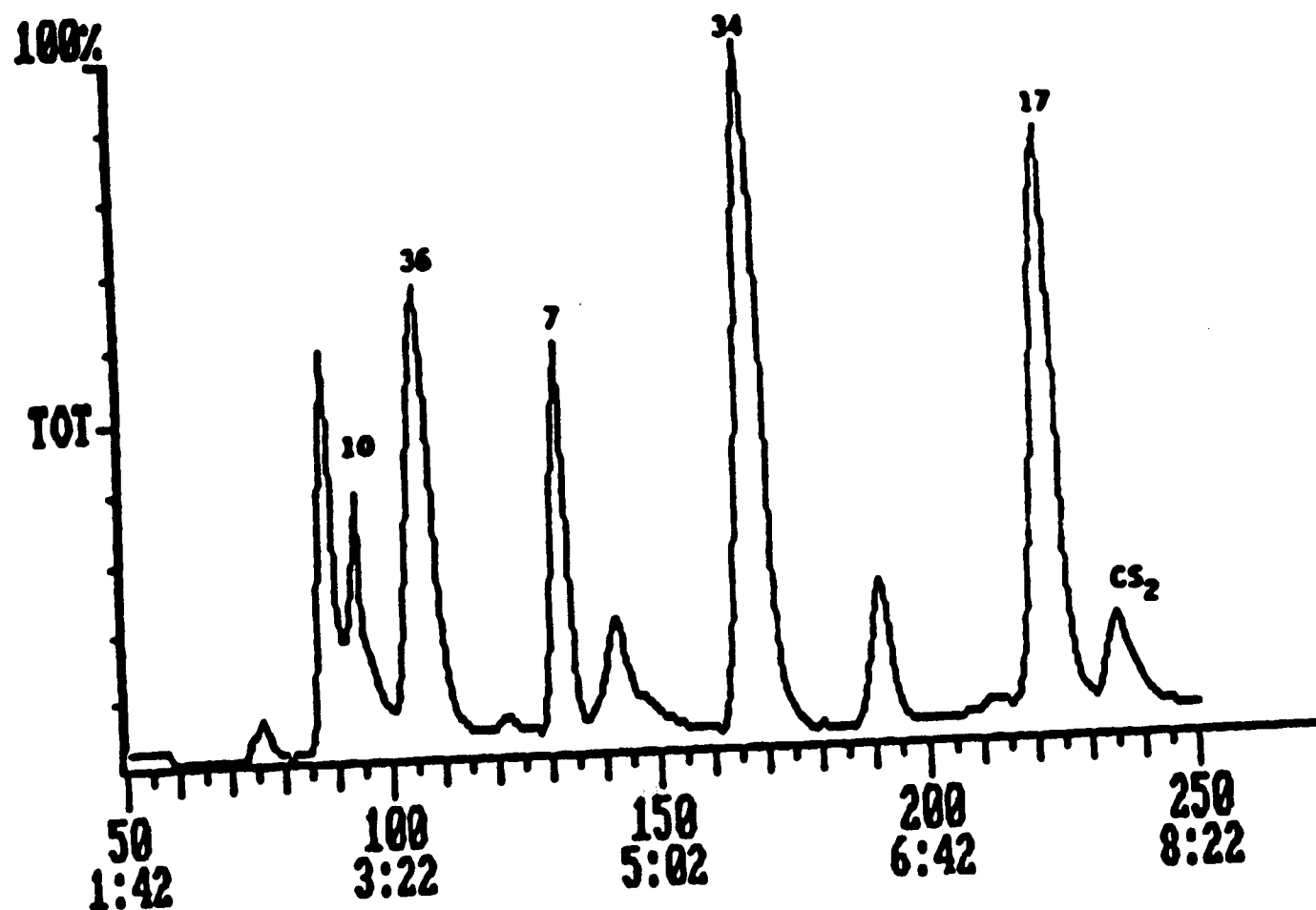



900  
2/25/91

**METHOD 200.8**

**DETERMINATION OF TRACE ELEMENTS IN WATERS AND WASTES  
BY INDUCTIVELY COUPLED PLASMA - MASS SPECTROMETRY**

**Stephen E. Long  
Technology Applications, Inc.**

**and**

**Theodore D. Martin  
Inorganic Chemistry Branch  
Chemistry Research Division**

**Revision 4.4  
April 1991**

**ENVIRONMENTAL MONITORING SYSTEMS LABORATORY  
OFFICE OF RESEARCH AND DEVELOPMENT  
U.S. ENVIRONMENTAL PROTECTION AGENCY  
CINCINNATI, OHIO 45268**

## METHOD 200.8

### DETERMINATION OF TRACE ELEMENTS IN WATERS AND WASTES BY INDUCTIVELY COUPLED PLASMA - MASS SPECTROMETRY

#### 1. SCOPE AND APPLICATION

- 1.1 This method provides procedures for determination of dissolved elements in ground waters, surface waters and drinking water. It may also be used for determination of total recoverable element concentrations in these waters as well as wastewaters, sludges and solid waste samples.
- 1.2 Dissolved elements are determined after suitable filtration and acid preservation. Acid digestion procedures are required prior to determination of total recoverable elements. In order to reduce potential interferences, dissolved solids should not exceed 0.2% (w/v) (Sect. 4.1.4).
- 1.3 This method is applicable to the following elements:

Element	Chemical Abstract Services Registry Numbers (CASRN)
Aluminum (Al)	7429-90-5
Antimony (Sb)	7440-36-0
Arsenic (As)	7440-38-2
Barium (Ba)	7440-39-3
Beryllium (Be)	7440-41-7
Cadmium (Cd)	7440-43-9
Chromium (Cr)	7440-47-3
Cobalt (Co)	7440-48-4
Copper (Cu)	7440-50-8
Lead (Pb)	7439-92-1
Manganese (Mn)	7439-96-5
Molybdenum (Mo)	7439-98-7
Nickel (Ni)	7440-02-0
Selenium (Se)	7782-49-2
Silver (Ag)	7440-22-4
Thallium (Tl)	7440-28-0
Thorium (Th)	7440-29-1
Uranium (U)	7440-61-1
Vanadium (V)	7440-62-2
Zinc (Zn)	7440-66-6

Estimated instrument detection limits (IDLs) for these elements are listed in Table 1. These are intended as a guide to instrumental limits typical of a system optimized for multielement determinations and employing commercial instrumentation and pneumatic nebulization sample introduction. However, actual method detection limits (MDLs)

and linear working ranges will be dependent on the sample matrix, instrumentation and selected operating conditions.

- 1.4 This method is suitable for the determination of silver in aqueous samples containing concentrations up to 0.1 mg/L. For the analysis of wastewater samples containing higher concentrations of silver, succeeding smaller volume, well mixed sample aliquots must be prepared until the analysis solution contains < 0.1 mg/L silver.
- 1.5 This method should be used by analysts experienced in the use of inductively coupled plasma mass spectrometry (ICP-MS), the interpretation of spectral and matrix interferences and procedures for their correction. A minimum of six months experience with commercial instrumentation is recommended.

## 2. SUMMARY OF METHOD

- 2.1 The method describes the multi-element determination of trace elements by ICP-MS<sup>1-3</sup>. Sample material in solution is introduced by pneumatic nebulization into a radiofrequency plasma where energy transfer processes cause desolvation, atomization and ionization. The ions are extracted from the plasma through a differentially pumped vacuum interface and separated on the basis of their mass-to-charge ratio by a quadrupole mass spectrometer having a minimum resolution capability of 1 amu peak width at 5% peak height. The ions transmitted through the quadrupole are registered by a continuous dynode electron multiplier or Faraday detector and the ion information processed by a data handling system. Interferences relating to the technique (Sect. 4) must be recognized and corrected for. Such corrections must include compensation for isobaric elemental interferences and interferences from polyatomic ions derived from the plasma gas, reagents or sample matrix. Instrumental drift as well as suppressions or enhancements of instrument response caused by the sample matrix must be corrected for by the use of internal standardization.

## 3. DEFINITIONS

- 3.1 DISSOLVED - Material that will pass through a 0.45  $\mu$ m membrane filter assembly, prior to sample acidification.
- 3.2 TOTAL RECOVERABLE - The concentration of analyte determined on an unfiltered sample following treatment with hot dilute mineral acid.
- 3.3 INSTRUMENT DETECTION LIMIT (IDL) - The concentration equivalent of the analyte signal, which is equal to three times the standard deviation of the blank signal at the selected analytical mass(es).
- 3.4 METHOD DETECTION LIMIT (MDL) - The minimum concentration of an analyte that can be identified, measured and reported with 99% confidence that the analyte concentration is greater than zero.

- 3.5 **LINEAR DYNAMIC RANGE (LDR)** - The concentration range over which the analytical working curve remains linear.
- 3.6 **LABORATORY REAGENT BLANK (LRB)** (preparation blank) - An aliquot of reagent water that is treated exactly as a sample including exposure to all labware, equipment, solvents, reagents, and internal standards that are used with other samples. The LRB is used to determine if method analytes or other interferences are present in the laboratory environment, the reagents or apparatus.
- 3.7 **CALIBRATION BLANK** - A volume of ASTM type I water acidified with the same acid matrix as is present in the calibration standards.
- 3.8 **INTERNAL STANDARD** - Pure analyte(s) added to a solution in known amount(s) and used to measure the relative responses of other method analytes that are components of the same solution. The internal standard must be an analyte that is not a sample component.
- 3.9 **STOCK STANDARD SOLUTION** - A concentrated solution containing one or more analytes prepared in the laboratory using assayed reference compounds or purchased from a reputable commercial source.
- 3.10 **CALIBRATION STANDARD (CAL)** - A solution prepared from the stock standard solution(s) which is used to calibrate the instrument response with respect to analyte concentration.
- 3.11 **TUNING SOLUTION** - A solution which is used to determine acceptable instrument performance prior to calibration and sample analyses.
- 3.12 **LABORATORY FORTIFIED BLANK (LFB)** - An aliquot of reagent water to which known quantities of the method analytes are added in the laboratory. The LFB is analyzed exactly like a sample, and its purpose is to determine whether method performance is within accepted control limits.
- 3.13 **LABORATORY FORTIFIED SAMPLE MATRIX (LFM)** - An aliquot of an environmental sample to which known quantities of the method analytes are added in the laboratory. The LFM is analyzed exactly like a sample, and its purpose is to determine whether the sample matrix contributes bias to the analytical results. The background concentrations of the analytes in the sample matrix must be determined in a separate aliquot and the measured values in the LFM corrected for the concentrations found.
- 3.14 **QUALITY CONTROL SAMPLE (QCS)** - A solution containing known concentrations of method analytes which is used to fortify an aliquot of LRB matrix. The QCS is obtained from a source external to the laboratory and is used to check laboratory performance.

#### **4. INTERFERENCES**

**4.1 Several interference sources may cause inaccuracies in the determination of trace elements by ICP-MS. These are:**

**4.1.1 Isobaric elemental interferences -** Are caused by isotopes of different elements which form singly or doubly charged ions of the same nominal mass-to-charge ratio and which cannot be resolved by the mass spectrometer in use. All elements determined by this method have, at a minimum, one isotope free of isobaric elemental interference. Of the analytical isotopes recommended for use with this method (Table 4), only molybdenum-98 (ruthenium) and selenium-82 (krypton) have isobaric elemental interferences. If alternative analytical isotopes having higher natural abundance are selected in order to achieve greater sensitivity, an isobaric interference may occur. All data obtained under such conditions must be corrected by measuring the signal from another isotope of the interfering element and subtracting the appropriate signal ratio from the isotope of interest. A record of this correction process should be included with the report of the data. It should be noted that such corrections will only be as accurate as the accuracy of the isotope ratio used in the elemental equation for data calculations. Relevant isotope ratios and instrument bias factors should be established prior to the application of any corrections.

**4.1.2 Abundance sensitivity -** Is a property defining the degree to which the wings of a mass peak contribute to adjacent masses. The abundance sensitivity is affected by ion energy and quadrupole operating pressure. Wing overlap interferences may result when a small ion peak is being measured adjacent to a large one. The potential for these interferences should be recognized and the spectrometer resolution adjusted to minimize them.

**4.1.3 Isobaric polyatomic ion interferences -** Are caused by ions consisting of more than one atom which have the same nominal mass-to-charge ratio as the isotope of interest, and which cannot be resolved by the mass spectrometer in use. These ions are commonly formed in the plasma or interface system from support gases or sample components. Most of the common interferences have been identified<sup>3</sup>, and these are listed in Table 2 together with the method elements affected. Such interferences must be recognized, and when they cannot be avoided by the selection of alternative analytical isotopes, appropriate corrections must be made to the data. Equations for the correction of data should be established at the time of the analytical run sequence as the polyatomic ion interferences will be highly dependent on the sample matrix and chosen instrument conditions.

4.1.4 Physical interferences - Are associated with the physical processes which govern the transport of sample into the plasma, sample conversion processes in the plasma, and the transmission of ions through the plasma-mass spectrometer interface. These interferences may result in differences between instrument responses for the sample and the calibration standards. Physical interferences may occur in the transfer of solution to the nebulizer (e.g., viscosity effects), at the point of aerosol formation and transport to the plasma (e.g., surface tension), or during excitation and ionization processes within the plasma itself. High levels of dissolved solids in the sample may contribute deposits of material on the extraction and/or skimmer cones reducing the effective diameter of the orifices and therefore ion transmission. Dissolved solids levels not exceeding 0.2% (w/v) have been recommended to reduce such effects. Internal standardization may be effectively used to compensate for many physical interference effects. Internal standards ideally should have similar analytical behavior to the elements being determined.

4.1.5 Memory interferences - Result when isotopes of elements in a previous sample contribute to the signals measured in a new sample. Memory effects can result from sample deposition on the sampler and skimmer cones, and from the buildup of sample material in the plasma torch and spray chamber. The site where these effects occur is dependent on the element and can be minimized by flushing the system with a rinse blank between samples (Sect. 7.6.3). The possibility of memory interferences should be recognized within an analytical run and suitable rinse times should be used to reduce them. The rinse times necessary for a particular element should be estimated prior to analysis. This may be achieved by aspirating a standard containing elements corresponding to ten times the upper end of the linear range for a normal sample analysis period, followed by analysis of the rinse blank at designated intervals. The length of time required to reduce analyte signals to within a factor of ten of the method detection limit, should be noted. Memory interferences may also be assessed within an analytical run by using a minimum of three replicate integrations for data acquisition. If the integrated signal values drop consecutively, the analyst should be alerted to the possibility of a memory effect, and should examine the analyte concentration in the previous sample to identify if this was high. If a memory interference is suspected, the sample should be reanalyzed after a long rinse period.

## 5. SAFETY

5.1 The toxicity or carcinogenicity of reagents used in this method have not been fully established. Each chemical should be regarded as a

potential health hazard and exposure to these compounds should be as low as reasonably achievable. Each laboratory is responsible for maintaining a current awareness file of OSHA regulations regarding the safe handling of the chemicals specified in this method<sup>6</sup>. A reference file of material data handling sheets should also be available to all personnel involved in the chemical analysis.

- 5.2 Analytical plasma sources emit radiofrequency radiation in addition to intense UV radiation. Suitable precautions should be taken to protect personnel from such hazards.

## **6. APPARATUS AND EQUIPMENT**

### **6.1 INDUCTIVELY COUPLED PLASMA - MASS SPECTROMETER**

- 6.1.1 Instrument capable of scanning the mass range 5-250 amu with a minimum resolution capability of 1 amu peak width at 5% peak height. Instrument may be fitted with a conventional or extended dynamic range detection system.
- 6.1.2 Argon gas supply (high-purity grade, 99.99%).
- 6.1.3 A variable-speed peristaltic pump is required for solution delivery to the nebulizer.
- 6.1.4 A mass-flow controller on the nebulizer gas supply is required. A water-cooled spray chamber may be of benefit in reducing some types of interferences (e.g., from polyatomic oxide species).
- 6.1.5 Operating conditions - Because of the diversity of instrument hardware, no detailed instrument operating conditions are provided. The analyst is advised to follow the recommended operating conditions provided by the manufacturer. It is the responsibility of the analyst to verify that the instrument configuration and operating conditions satisfy the analytical requirements and to maintain quality control data verifying instrument performance and analytical results. Instrument operating conditions which were used to generate precision and recovery data for this method (Sect. 13) are included in Table 6.
- 6.1.6 If an electron multiplier detector is being used, precautions should be taken, where necessary, to prevent exposure to high ion flux. Otherwise changes in instrument response or damage to the multiplier may result. Samples having high concentrations of elements beyond the linear range of the instrument and with isotopes falling within scanning windows should be diluted prior to analysis.
- 6.2 LABWARE - For the determination of trace levels of elements, contamination and loss are of prime consideration. Potential

contamination sources include improperly cleaned laboratory apparatus and general contamination within the laboratory environment from dust, etc. A clean laboratory work area, designated for trace element sample handling must be used. Sample containers can introduce positive and negative errors in the determination of trace elements by (1) contributing contaminants through surface desorption or leaching, (2) depleting element concentrations through adsorption processes. All reuseable labware (glass, quartz, polyethylene, Teflon, etc.) including the sample container should be cleaned prior to use. Labware may be soaked overnight and thoroughly washed with laboratory-grade detergent and water, rinsed with water, and soaked for four hours in a mixture of dilute nitric and hydrochloric acid (1+2+9), followed by rinsing with water, ASTM type I water and oven drying.

**NOTE:** Chromic acid must not be used for cleaning glassware.

6.2.1 Glassware - Volumetric flasks, graduated cylinders, funnels and centrifuge tubes.

6.2.2 Assorted calibrated pipettes.

6.2.3 Conical Phillips beakers, 250-mL with 50-mm watch glasses.  
Griffin beakers, 250-mL with 75-mm watch glasses.

6.2.4 Storage bottles - Narrow mouth bottles, Teflon FEP (fluorinated ethylene propylene) with Tefzel ETFE (ethylene tetrafluorethylene) screw closure, 125-mL and 250-mL capacities.

### 6.3 SAMPLE PROCESSING EQUIPMENT

6.3.1 Air Displacement Pipetter - Digital pipet system capable of delivering volumes from 10 to 2500  $\mu$ L with an assortment of high quality disposable pipet tips.

6.3.2 Balance - Analytical, capable of accurately weighing to 0.1 mg.

6.3.3 Hot Plate - (Corning PC100 or equivalent).

6.3.4 Centrifuge - Steel cabinet with guard bowl, electric timer and brake.

6.3.5 Drying Oven - Gravity convection oven with thermostatic control capable of maintaining  $105^{\circ}\text{C} \pm 5^{\circ}\text{C}$ .

## 7. REAGENTS AND CONSUMABLE MATERIALS

7.1 Reagents may contain elemental impurities that might affect the integrity of analytical data. Owing to the high sensitivity of ICP-



MS, high-purity reagents should be used whenever possible. All acids used for this method must be of ultra high-purity grade. Suitable acids are available from a number of manufacturers or may be prepared by sub-boiling distillation. Nitric acid is preferred for ICP-MS in order to minimize polyatomic ion interferences. Several polyatomic ion interferences result when hydrochloric acid is used (Table 2), however, it should be noted that hydrochloric acid is required to maintain stability in solutions containing antimony and silver. When hydrochloric acid is used, corrections for the chloride polyatomic ion interferences must be applied to all data.

- 7.1.1 Nitric acid, concentrated (sp.gr. 1.41).
- 7.1.2 Nitric acid (1+1) - Add 500 mL conc. nitric acid to 400 mL of ASTM type I water and dilute to 1 L.
- 7.1.3 Nitric acid (1+9) - Add 100 mL conc. nitric acid to 400 mL of ASTM type I water and dilute to 1 L.
- 7.1.4 Hydrochloric acid, concentrated (sp.gr. 1.19).
- 7.1.5 Hydrochloric acid (1+1) - Add 500 mL conc. hydrochloric acid to 400 mL of ASTM type I water and dilute to 1 L.
- 7.1.6 Hydrochloric acid (1+4) - Add 200 mL conc. hydrochloric acid to 400 mL of ASTM type I water and dilute to 1 L.
- 7.1.7 Ammonium hydroxide, concentrated (sp.gr. 0.902).
- 7.1.8 Tartaric acid (CASRN 87-69-4).
- 7.2 WATER - For all sample preparation and dilutions, ASTM type I water (ASTM D1193) is required. Suitable water may be prepared by passing distilled water through a mixed bed of anion and cation exchange resins.
- 7.3 STANDARD STOCK SOLUTIONS - May be purchased from a reputable commercial source or prepared from ultra high-purity grade chemicals or metals (99.99 - 99.999% pure). All salts should be dried for 1 h at 105°C, unless otherwise specified. (CAUTION: Many metal salts are extremely toxic if inhaled or swallowed. Wash hands thoroughly after handling). Stock solutions should be stored in Teflon bottles. The following procedures may be used for preparing standard stock solutions:

NOTE: Some metals, particularly those which form surface oxides require cleaning prior to being weighed. This may be achieved by pickling the surface of the metal in acid. An amount in excess of the desired weight should be pickled repeatedly, rinsed with water, dried and weighed until the desired weight is achieved.

- 7.3.1 Aluminum solution, stock 1 mL = 1000  $\mu$ g Al: Pickle aluminum metal in warm (1+1) HCl to an exact weight of 0.100 g. Dissolve in 10 mL conc. HCl and 2 mL conc. nitric acid, heating to effect solution. Continue heating until volume is reduced to 4 mL. Cool and add 4 mL ASTM type I water. Heat until the volume is reduced to 2 mL. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.2 Antimony solution, stock 1 mL = 1000  $\mu$ g Sb: Dissolve 0.100 g antimony powder in 2 mL (1+1) nitric acid and 0.5 mL conc. hydrochloric acid, heating to effect solution. Cool, add 20 mL ASTM type I water and 0.15 g tartaric acid. Warm the solution to dissolve the white precipitate. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.3 Arsenic solution, stock 1 mL = 1000  $\mu$ g As: Dissolve 0.1320 g  $\text{As}_2\text{O}_3$  in a mixture of 50 mL ASTM type I water and 1 mL conc. ammonium hydroxide. Heat gently to dissolve. Cool and acidify the solution with 2 mL conc. nitric acid. Dilute to 100 mL with ASTM type I water.
- 7.3.4 Barium solution, stock 1 mL = 1000  $\mu$ g Ba: Dissolve 0.1437 g  $\text{BaCO}_3$  in a solution mixture of 10 mL ASTM type I water and 2 mL conc. nitric acid. Heat and stir to effect solution and degassing. Dilute to 100 mL with ASTM type I water.
- 7.3.5 Beryllium solution, stock 1 mL = 1000  $\mu$ g Be: Dissolve 1.965 g  $\text{BeSO}_4 \cdot 4\text{H}_2\text{O}$  (DO NOT DRY) in 50 mL ASTM Type I water. Add 1 mL conc. nitric acid. Dilute to 100 mL with ASTM type I water.
- 7.3.6 Bismuth solution, stock 1 mL = 1000  $\mu$ g Bi: Dissolve 0.1115 g  $\text{Bi}_2\text{O}_3$  in 5 mL conc. nitric acid. Heat to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.7 Cadmium solution, stock 1 mL = 1000  $\mu$ g Cd: Pickle cadmium metal in (1+9) nitric acid to an exact weight of 0.100 g. Dissolve in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.8 Chromium solution, stock 1 mL = 1000  $\mu$ g Cr: Dissolve 0.1923 g  $\text{CrO}_3$  in a solution mixture of 10 mL ASTM type I water and 1 mL conc. nitric acid. Dilute to 100 mL with ASTM type I water.
- 7.3.9 Cobalt solution, stock 1 mL = 1000  $\mu$ g Co: Pickle cobalt metal in (1+9) nitric acid to an exact weight of 0.100 g. Dissolve in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.

- 7.3.10 Copper solution, stock 1 mL = 1000  $\mu$ g Cu: Pickle copper metal in (1+9) nitric acid to an exact weight of 0.100 g. Dissolve in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.11 Indium solution, stock 1 mL = 1000  $\mu$ g In: Pickle indium metal in (1+1) nitric acid to an exact weight of 0.100 g. Dissolve in 10 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.12 Lead solution, stock 1 mL = 1000  $\mu$ g Pb: Dissolve 0.1599 g  $\text{PbNO}_3$  in 5 mL (1+1) nitric acid. Dilute to 100 mL with ASTM type I water.
- 7.3.13 Magnesium solution, stock 1 mL = 1000  $\mu$ g Mg: Dissolve 0.1658 g  $\text{MgO}$  in 10 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.14 Manganese solution, stock 1 mL = 1000  $\mu$ g Mn: Pickle manganese flake in (1+9) nitric acid to an exact weight of 0.100 g. Dissolve in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.15 Molybdenum solution, stock 1 mL = 1000  $\mu$ g Mo: Dissolve 0.1500 g  $\text{MoO}_3$  in a solution mixture of 10 mL ASTM type I water and 1 mL conc. ammonium hydroxide., heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.16 Nickel solution, stock 1 mL = 1000  $\mu$ g Ni: Dissolve 0.100 g nickel powder in 5 mL conc. nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.17 Scandium solution, stock 1 mL = 1000  $\mu$ g Sc: Dissolve 0.1534 g  $\text{Sc}_2\text{O}_3$  in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.18 Selenium solution, stock 1 mL = 1000  $\mu$ g Se: Dissolve 0.1405 g  $\text{SeO}_2$  in 20 mL ASTM type I water. Dilute to 100 mL with ASTM type I water.
- 7.3.19 Silver solution, stock 1 mL = 1000  $\mu$ g Ag: Dissolve 0.100 g silver metal in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water. Store in dark container.
- 7.3.20 Terbium solution, stock 1 mL = 1000  $\mu$ g Tb: Dissolve 0.1176 g  $\text{Tb}_4\text{O}_7$  in 5 mL conc. nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.

- 7.3.21 Thallium solution, stock 1 mL = 1000  $\mu$ g Tl: Dissolve 0.1303 g  $\text{TlNO}_3$  in a solution mixture of 10 mL ASTM type I water and 1 mL conc. nitric acid. Dilute to 100 mL with ASTM type I water.
- 7.3.22 Thorium solution, stock 1 mL = 1000  $\mu$ g Th: Dissolve 0.2380 g  $\text{Th}(\text{NO}_3)_4 \cdot 4\text{H}_2\text{O}$  (DO NOT DRY) in 20 mL ASTM type I water. Dilute to 100 mL with ASTM type I water.
- 7.3.23 Uranium solution, stock 1 mL = 1000  $\mu$ g U: Dissolve 0.2110 g  $\text{UO}_2(\text{NO}_3)_2 \cdot 6\text{H}_2\text{O}$  (DO NOT DRY) in 20 mL ASTM type I water and dilute to 100 mL with ASTM type I water.
- 7.3.24 Vanadium solution, stock 1 mL = 1000  $\mu$ g V: Pickle vanadium metal in (1+9) nitric acid to an exact weight of 0.100 g. Dissolve in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.25 Yttrium solution, stock 1 mL = 1000  $\mu$ g Y: Dissolve 0.1270 g  $\text{Y}_2\text{O}_3$  in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.3.26 Zinc solution, stock 1 mL = 1000  $\mu$ g Zn: Pickle zinc metal in (1+9) nitric acid to an exact weight of 0.100 g. Dissolve in 5 mL (1+1) nitric acid, heating to effect solution. Cool and dilute to 100 mL with ASTM type I water.
- 7.4 MULTIELEMENT STOCK STANDARD SOLUTIONS - Care must be taken in the preparation of multielement stock standards that the elements are compatible and stable. Originating element stocks should be checked for the presence of impurities which might influence the accuracy of the standard. Freshly prepared standards should be transferred to acid cleaned, not previously used FEP fluorocarbon bottles for storage and monitored periodically for stability. The following combinations of elements are suggested:

#### Standard Solution A

Aluminum	Manganese
Antimony	Molybdenum
Arsenic	Nickel
Beryllium	Selenium
Cadmium	Thallium
Chromium	Thorium
Cobalt	Uranium
Copper	Vanadium
Lead	Zinc

#### Standard Solution B

Barium
Silver

Multielement stock standard solutions A and B (1 mL = 10  $\mu$ g) may be prepared by diluting 1 mL of each single element stock in the

combination list to 100 mL with ASTM type I water containing 1% (v/v) nitric acid.

7.4.1 Preparation of calibration standards - fresh multielement calibration standards should be prepared every two weeks or as needed. Dilute each of the stock multielement standard solutions A and B to levels appropriate to the operating range of the instrument using ASTM type I water containing 1% (v/v) nitric acid. The element concentrations in the standards should be sufficiently high to produce good measurement precision and to accurately define the slope of the response curve. Concentrations of 200 µg/L are suggested. If the direct addition procedure is being used (Method A, Sect. 9.2), add internal standards (Sect. 7.5) to the calibration standards and store in Teflon bottles. Calibration standards should be verified initially using a quality control sample (Sect. 7.8).

7.5 INTERNAL STANDARDS STOCK SOLUTION, 1 mL = 100 µg. Dilute 10 mL of scandium, yttrium, indium, terbium and bismuth stock standards (Sect. 7.3) to 100 mL with ASTM type I water, and store in a Teflon bottle. Use this solution concentrate for addition to blanks, calibration standards and samples, or dilute by an appropriate amount using 1% (v/v) nitric acid, if the internal standards are being added by peristaltic pump (Method B, Sect. 9.2).

7.6 BLANKS - Three types of blanks are required for this method. A calibration blank is used to establish the analytical calibration curve, the laboratory reagent blank is used to assess possible contamination from the sample preparation procedure and to assess spectral background and the rinse blank is used to flush the instrument between samples in order to reduce memory interferences.

7.6.1 Calibration blank - Consists of 1% (v/v) nitric acid in ASTM type I water. If the direct addition procedure (Method A, Sect. 9.2), is being used add internal standards.

7.6.2 Laboratory reagent blank (LRB) - Must contain all the reagents in the same volumes as used in processing the samples. The LRB must be carried through the entire sample digestion and preparation scheme. If the direct addition procedure (Method A, Sect. 9.2) is being used, add internal standards to the solution after preparation is complete.

7.6.3 Rinse blank - Consists of 2% (v/v) nitric acid in ASTM type I water.

7.7 TUNING SOLUTION - This solution is used for instrument tuning and mass calibration prior to analysis. The solution is prepared by mixing beryllium, magnesium, cobalt, indium and lead stock solutions (Sect. 7.3) in 1% (v/v) nitric acid to produce a concentration of

100 µg/L of each element. Internal standards are not added to this solution.

7.8 **QUALITY CONTROL SAMPLE (QCS)** - The QCS should be obtained from a source outside the laboratory. Dilute an appropriate aliquot of analytes (concentrations not to exceed 1000 µg/L), in 1% (v/v) nitric acid. If the direct addition procedure (Method A, Sect. 9.2) is being used, add internal standards after dilution, mix and store in a Teflon bottle.

7.9 **LABORATORY FORTIFIED BLANK (LFB)** - To an aliquot of LRB, add aliquots from multielement stock standards A and B (Sect. 7.4) to produce a final concentration of 100 µg/L for each analyte. The LFB must be carried through the entire sample digestion and preparation scheme. If the direct addition procedure (Method A, Sect. 9.2) is being used, add internal standards to this solution after preparation has been completed.

## **8. SAMPLE COLLECTION, PRESERVATION AND STORAGE**

8.1 Prior to sample collection, consideration should be given to the type of data required so that appropriate preservation and pretreatment steps can be taken. Filtration, acid preservation, etc., should be performed at the time of sample collection or as soon thereafter as practically possible.

8.2 For the determination of dissolved elements, the sample should be filtered through a 0.45-µm membrane filter. Use a portion of the sample to rinse the filter assembly, discard and then collect the required volume of filtrate. Acidify the filtrate with (1+1) nitric acid immediately following filtration to pH < 2.

8.3 For the determination of total recoverable elements in aqueous samples, acidify with (1+1) nitric acid at the time of collection to pH < 2 (normally, 3 mL of (1+1) nitric acid per liter of sample is sufficient for most ambient and drinking water samples). The sample should not be filtered prior to analysis.

**NOTE:** Samples that cannot be acid preserved at the time of collection because of sampling limitations or transport restrictions, should be acidified with nitric acid to pH < 2 upon receipt in the laboratory. Following acidification, the sample should be held for 16 h before withdrawing an aliquot for sample processing.

8.4 Solid samples usually require no preservation prior to analysis other than storage at 4°C.

## **9. CALIBRATION AND STANDARDIZATION**

**9.1 CALIBRATION** - Demonstration and documentation of acceptable initial calibration is required before any samples are analyzed and is required periodically throughout sample analysis as dictated by results of continuing calibration checks. After initial calibration is successful, a calibration check is required at the beginning and end of each period during which analyses are performed, and at requisite intervals.

**9.1.1** Initiate proper operating configuration of instrument and data system. Allow a period of not less than 30 min for the instrument to warm up. During this process conduct mass calibration and resolution checks using the tuning solution. Resolution at low mass is indicated by magnesium isotopes 24,25,26. Resolution at high mass is indicated by lead isotopes 206,207,208. For good performance adjust spectrometer resolution to produce a peak width of approximately 0.75 amu at 5% peak height. Adjust mass calibration if it has shifted by more than 0.1 amu from unit mass.

**9.1.2** Instrument stability must be demonstrated by running the tuning solution (Sect. 7.7) a minimum of five times with resulting relative standard deviations of absolute signals for all analytes of less than 5%.

**9.1.3** Prior to initial calibration, set up proper instrument software routines for quantitative analysis. The instrument must be calibrated for the analytes to be determined using the calibration blank (Sect. 7.6.1) and calibration standards A and B (Sect. 7.4.1) prepared at one or more concentration levels. A minimum of three replicate integrations are required for data acquisition. Use the average of the integrations for instrument calibration and data reporting.

**9.1.4** The rinse blank should be used to flush the system between solution changes for blanks, standards and samples. Allow sufficient rinse time to remove traces of the previous sample or a minimum of 1 min. Solutions should be aspirated for 30 sec prior to the acquisition of data to allow equilibrium to be established.

**9.2 INTERNAL STANDARDIZATION** - Internal standardization must be used in all analyses to correct for instrument drift and physical interferences. A list of acceptable internal standards is provided in Table 3. For full mass range scans, a minimum of three internal standards must be used. Procedures described in this method for general application, detail the use of five internal standards; scandium, yttrium, indium, terbium and bismuth. These were used to generate the precision and recovery data attached to this method. Internal standards must be present in all samples, standards and

blanks at identical levels. This may be achieved by directly adding an aliquot of the internal standards to the CAL standard, blank or sample solution (Method A, Sect. 9.2), or alternatively by mixing with the solution prior to nebulization using a second channel of the peristaltic pump and a mixing coil (Method B, Sect. 9.2). The concentration of the internal standard should be sufficiently high that good precision is obtained in the measurement of the isotope used for data correction and to minimize the possibility of correction errors if the internal standard is naturally present in the sample. A concentration of 200 µg/L of each internal standard is recommended. Internal standards should be added to blanks, samples and standards in a like manner, so that dilution effects resulting from the addition may be disregarded.

**9.3 INSTRUMENT PERFORMANCE** - Check the performance of the instrument and verify the calibration using data gathered from analyses of calibration blanks, calibration standards and the quality control sample (QCS).

9.3.1 After the calibration has been established, it must be initially verified for all analytes by analyzing the QCS (Sect. 7.8). If measurements exceed  $\pm 10\%$  of the established QCS value, the analysis should be terminated, the source of the problem identified and corrected, the instrument recalibrated and the calibration reverified before continuing analyses.

9.3.2 To verify that the instrument is properly calibrated on a continuing basis, run the calibration blank and calibration standards as surrogate samples after every ten analyses. The results of the analyses of the standards will indicate whether the calibration remains valid. If the indicated concentration of any analyte deviates from the true concentration by more than 10%, reanalyze the standard. If the analyte is again outside the 10% limit, the instrument must be recalibrated and the previous ten samples reanalyzed. The instrument responses from the calibration check may be used for recalibration purposes. If the sample matrix is responsible for the calibration drift, it is recommended that the previous ten samples are reanalyzed in groups of five between calibration checks to prevent a similar drift situation from occurring.

**10. QUALITY CONTROL**

10.1 Each laboratory using this method is required to operate a formal quality control (QC) program. The minimum requirements of this program consist of an initial demonstration of laboratory capability, and the analysis of laboratory reagent blanks, fortified blanks and samples as a continuing check on performance. The



laboratory is required to maintain performance records that define the quality of the data thus generated.

## 10.2 INITIAL DEMONSTRATION OF PERFORMANCE

10.2.1 The initial demonstration of performance is used to characterize instrument performance (method detection limits and linear calibration ranges) for analyses conducted by this method.

10.2.2 Method detection limits (MDL) should be established for all analytes, using reagent water (blank) fortified at a concentration of two to five times the estimated detection limit. To determine MDL values, take seven replicate aliquots of the fortified reagent water and process through the entire analytical method. Perform all calculations defined in the method and report the concentration values in the appropriate units. Calculate the MDL as follows:

$$MDL = (t) \times (S)$$

where,  $t$  = Student's  $t$  value for a 99% confidence level and a standard deviation estimate with  $n-1$  degrees of freedom [ $t = 3.14$  for seven replicates].

$S$  = standard deviation of the replicate analyses.

MDLs should be determined every six months or whenever a significant change in background or instrument response is expected (e.g., detector change).

10.2.3 Linear calibration ranges - Linear calibration ranges are primarily detector limited. The upper limit of the linear calibration range should be established for each analyte by determining the signal responses from a minimum of three different concentration standards, one of which is close to the upper limit of the linear range. Care should be taken to avoid potential damage to the detector during this process. The linear calibration range which may be used for the analysis of samples should be judged by the analyst from the resulting data. Linear calibration ranges should be determined every six months or whenever a significant change in instrument response is expected (e.g., detector change).

## 10.3 ASSESSING LABORATORY PERFORMANCE - REAGENT AND FORTIFIED BLANKS

10.3.1 Laboratory reagent blank (LRB) - The laboratory must analyze at least one LRB (Sect. 7.6.2) with each set of samples. LRB data are used to assess contamination from the laboratory environment and to characterize spectral background from the reagents used in sample processing. If an analyte value in the reagent blank exceeds its determined MDL, then laboratory

or reagent contamination should be suspected. Any determined source of contamination should be corrected and the samples reanalyzed.

10.3.2 Laboratory fortified blank (LFB) - The laboratory must analyze at least one LFB (Sect. 7.9) with each batch of samples. Calculate accuracy as percent recovery (Sect. 10.4.2). If the recovery of any analyte falls outside the control limits (Sect. 10.3.3), that analyte is judged out of control, and the source of the problem should be identified and resolved before continuing analyses.

10.3.3 Until sufficient LFB data become available (usually a minimum of 20 to 30 analyses), the laboratory should assess laboratory performance against recovery limits of 85-115%. When sufficient internal performance data becomes available, develop control limits from the percent mean recovery ( $\bar{x}$ ) and the standard deviation ( $S$ ) of the mean recovery. These data are used to establish upper and lower control limits as follows:

$$\begin{aligned}\text{UPPER CONTROL LIMIT} &= \bar{x} + 3S \\ \text{LOWER CONTROL LIMIT} &= \bar{x} - 3S\end{aligned}$$

After each five to ten new recovery measurements, new control limits should be calculated using only the most recent twenty to 30 data points.

#### 10.4 ASSESSING ANALYTE RECOVERY - LABORATORY FORTIFIED SAMPLE MATRIX

10.4.1 The laboratory must add a known amount of analyte to a minimum of 10% of the routine samples or one sample per sample set, whichever is greater. Ideally for water samples, the analyte concentration should be the same as that used in the LFB (Sect. 10.3.2). For solid samples, the concentration added should be 50 mg/kg equivalent (100  $\mu\text{g/L}$  in the analysis solution). Over time, samples from all routine sample sources should be fortified.

10.4.2 Calculate the percent recovery for each analyte, corrected for background concentrations measured in the unfortified sample, and compare these values to the control limits established in Sect. 10.3.3 for the analyses of LFBs. Recovery calculations are not required if the concentration of the analyte added is less than 10% of the sample background concentration. Percent recovery may be calculated in units appropriate to the matrix, using the following equation:

$$R = \frac{C_s - C}{s} \times 100$$

where, R = percent recovery  
C = fortified sample concentration  
C<sup>s</sup> = sample background concentration  
s = concentration equivalent of  
fortifier added to sample.

10.4.3 If recovery of any analyte falls outside the designated range and laboratory performance for that analyte is shown to be in control (Sect. 10.3), the recovery problem encountered with the fortified sample is judged to be matrix related, not system related. The result for that analyte in the unfortified sample must be labelled "suspect/matrix" to inform the data user that the results are suspect due to matrix effects.

10.5 INTERNAL STANDARDS RESPONSES - The analyst is expected to monitor the responses from the internal standards throughout the sample set being analyzed. Ratios of the internal standards responses against each other should also be monitored routinely. This information may be used to detect potential problems caused by mass dependent drift, errors incurred in adding the internal standards or increases in the concentrations of individual internal standards caused by background contributions from the sample. The absolute response of any one internal standard should not deviate more than 60-125% of the original response in the calibration blank. If deviations greater than this are observed, use the following test procedure:

10.5.1 Flush the instrument with the rinse blank and monitor the responses in the calibration blank. If the responses of the internal standards are now within the limit, take a fresh aliquot of the sample, dilute by a further factor of two, add the internal standards and reanalyze.

10.5.2 If test (Sect. 10.5.1) is not satisfied, or if it is a blank or calibration standard that is out of limits, terminate the analysis, and determine the cause of the drift. Possible causes of drift may be a partially blocked sampling cone or a change in the tuning condition of the instrument.

## 11. PROCEDURE

### 11.1 SAMPLE PREPARATION - DISSOLVED ELEMENTS

11.1.1 For determination of dissolved elements in drinking water, ground and surface waters, take a 100 mL aliquot of the filtered acid preserved sample, and add 1 mL of concentrated nitric acid. If the direct addition procedure (Method A) is being used, add internal standards and mix. The sample is now ready for analysis. Allowance for sample dilution should be made in the calculations.

**NOTE:** If a precipitate is formed during acidification, transport or storage, the sample aliquot must be treated using the procedure in Sect. 11.2.1 prior to analysis.

## **11.2 SAMPLE PREPARATION - TOTAL RECOVERABLE ELEMENTS**

**11.2.1** For determination of total recoverable elements in water or wastewater, take a 100 mL aliquot from a well mixed, acid preserved sample containing not more than 0.25% (w/v) total solids and transfer to a 250-mL Griffin beaker (if total solids are greater than 0.25% reduce the size of the aliquot by a proportionate amount). Add 1 mL of conc. nitric acid and 0.5 mL conc. hydrochloric acid. Heat on a hot plate at 85°C until the volume has been reduced to approximately 20 mL, ensuring that the sample does not boil. A spare beaker containing 20 mL of water can be used as a gauge. (NOTE: Adjust the temperature control of the hot plate such that an uncovered beaker containing 50 mL of water located in the center of the hot plate can be maintained at a temperature no higher than 85°C. Evaporation time for 100 mL of sample at 85°C is approximately 2 h with the rate of evaporation increasing rapidly as the sample volume approaches 20 mL). Cover the beaker with a watch glass and reflux for 30 min. Slight boiling may occur but vigorous boiling should be avoided. Allow to cool and quantitatively transfer to either a 50-mL volumetric flask or 50-mL class A stoppered graduated cylinder. Dilute to volume with ASTM type I water and mix. Centrifuge the sample or allow to stand overnight to separate insoluble material. Prior to analysis, pipette 20 mL into a 50-mL volumetric flask, dilute to volume with ASTM type I water and mix. If the direct addition procedure (Method A, Sect. 9.2) is being used, add internal standards and mix. The sample is now ready for analysis. Because the stability of diluted samples cannot be fully characterized, all analyses should be performed as soon as possible after the completed preparation.

**11.2.2** For determination of total recoverable elements in solid samples (sludge, soils, and sediments), mix the sample thoroughly to achieve homogeneity and weigh accurately a  $1.0 \pm 0.01$  g portion of the sample. Transfer to a 250-mL Phillips beaker. Add 4 mL (1+1) nitric acid and 10 mL (1+4) HCl. Cover with a watch glass, and reflux the sample on a hot plate for 30 min. Very slight boiling may occur, however, vigorous boiling must be avoided to prevent the loss of the HCl-H<sub>2</sub>O azeotrope. (NOTE: Adjust the temperature control of the hot plate such that an uncovered Griffin beaker containing 50 mL of water located in the center of the hot plate can be maintained at a temperature of approximately but no higher than 85°C). Allow the sample to cool, and quantitatively transfer to a 100-mL volumetric flask. Dilute

to volume with ASTM type I water and mix. Centrifuge the sample or allow to stand overnight to separate insoluble material. Prior to analysis, pipette 10 mL into a 50-mL volumetric flask and dilute to volume with ASTM type I water. If the direct addition procedure (Method A, Sect. 9.2) is being used, add internal standards and mix. The sample is now ready for analysis. Because the effects of various matrices on the stability of diluted samples cannot be characterized, all analyses should be performed as soon as possible after the completed preparation.

**NOTE:** Determine the percent solids in the sample for use in calculations and for reporting data on a dry weight basis.

- 11.3 For every new or unusual matrix, it is highly recommended that a semi-quantitative analysis be carried out to screen for high element concentrations. Information gained from this may be used to prevent potential damage to the detector during sample analysis and to identify elements which may be higher than the linear range. Matrix screening may be carried out by using intelligent software, if available, or by diluting the sample by a factor of 500 and analyzing in a semi-quantitative mode. The sample should also be screened for background levels of all elements chosen for use as internal standards in order to prevent bias in the calculation of the analytical data.
- 11.4 Initiate instrument operating configuration. Tune and calibrate the instrument for the analytes of interest (Sect. 9).
- 11.5 Establish instrument software run procedures for quantitative analysis. For all sample analyses, a minimum of three replicate integrations are required for data acquisition. Discard any integrations which are considered to be statistical outliers and use the average of the integrations for data reporting.
- 11.6 All masses which might affect data quality must be monitored during the analytical run. As a minimum, those masses prescribed in Table 4 must be monitored in the same scan as is used for the collection of the data. This information should be used to correct the data for identified interferences.
- 11.7 The rinse blank should be used to flush the system between samples. Allow sufficient time to remove traces of the previous sample or a minimum of one minute. Samples should be aspirated for 30 sec prior to the collection of data.
- 11.8 Samples having concentrations higher than the established linear dynamic range should be diluted into range and reanalyzed. The sample should first be analyzed for the trace elements in the sample, protecting the detector from the high concentration elements, if necessary, by the selection of appropriate scanning windows. The sample should then be diluted for the determination of

the remaining elements. Alternatively, the dynamic range may be adjusted by selecting an alternative isotope of lower natural abundance, provided quality control data for that isotope have been established. The dynamic range must not be adjusted by altering instrument conditions to an uncharacterized state.

## **12. CALCULATIONS**

- 12.1 Elemental equations recommended for sample data calculations are listed in Table 5. Sample data should be reported in units of  $\mu\text{g/L}$  for aqueous samples or  $\text{mg/kg}$  dry weight for solid samples. Do not report element concentrations below the determined MDL.
- 12.2 For data values less than ten, two significant figures should be used for reporting element concentrations. For data values greater than or equal to ten, three significant figures should be used.
- 12.3 Reported values should be calibration blank subtracted. For aqueous samples prepared by total recoverable procedure (Sect. 11.2.1), multiply solution concentrations by the dilution factor 1.25. For solid samples prepared by total recoverable procedure (Sect. 11.2.2), multiply solution concentrations ( $\mu\text{g/L}$  in the analysis solution) by the dilution factor 0.5. If additional dilutions were made to any samples, the appropriate factor should be applied to the calculated sample concentrations.
- 12.4 Data values should be corrected for instrument drift or sample matrix induced interferences by the application of internal standardization. Corrections for characterized spectral interferences should be applied to the data. Chloride interference corrections should be made on all samples, regardless of the addition of hydrochloric acid, as the chloride ion is a common constituent of environmental samples.
- 12.5 If an element has more than one monitored isotope, examination of the concentration calculated for each isotope, or the isotope ratios, will provide useful information for the analyst in detecting a possible spectral interference. Consideration should therefore be given to both primary and secondary isotopes in the evaluation of the element concentration. In some cases, secondary isotopes may be less sensitive or more prone to interferences than the primary recommended isotopes, therefore differences between the results do not necessarily indicate a problem with data calculated for the primary isotopes.
- 12.6 The QC data obtained during the analyses provide an indication of the quality of the sample data and should be provided with the sample results.

### **13. PRECISION AND ACCURACY**

- 13.1 Instrument operating conditions used for single laboratory testing of the method are summarized in Table 6. Total recoverable MDLs determined using the procedure described in Sect. 10.2.2, are listed in Table 7.
- 13.2 Data obtained from single laboratory testing of the method are summarized in Table 8 for five water samples representing drinking water, surface water, ground water and waste effluent. Samples were prepared using the procedure described in Sect. 11.2.1. For each matrix, five replicates were analyzed and the average of the replicates used for determining the sample background concentration for each element. Two further pairs of duplicates were fortified at different concentration levels. For each method element, the sample background concentration, mean percent recovery, the standard deviation of the percent recovery and the relative percent difference between the duplicate fortified samples are listed in Table 8.
- 13.3 Data obtained from single laboratory testing of the method are summarized in Table 9 for three solid samples consisting of SRM 1645 River Sediment, EPA Hazardous Soil and EPA Electroplating Sludge. Samples were prepared using the procedure described in Sect. 11.2.2. For each method element, the sample background concentration, mean percent recovery, the standard deviation of the percent recovery and the relative percent difference between the duplicate fortified samples were determined as for Sect. 13.2.

### **14. REFERENCES**

1. A. L. Gray and A. R. Date, *Analyst* 108 1033 (1983).
2. R. S. Houk et al. *Anal Chem.* 52 2283 (1980).
3. R. S. Houk, *Anal. Chem.* 58 97A (1986).
4. J. J. Thompson and R. S. Houk, *Appl. Spec.* 41 801 (1987).
5. "OSHA Safety and Health Standards, General Industry," (29 CFR 1910), Occupational Safety and Health Administration, OSHA 2206, revised January 1976.
6. "Proposed OSHA Safety and Health Standards, Laboratories," Occupational Safety and Health Administration, Federal Register, July 24, 1986.
7. Code of Federal Regulations 40, Ch. 1, Pt. 136 Appendix B.

**TABLE 1: ESTIMATED INSTRUMENT DETECTION LIMITS**

<b>ELEMENT</b>	<b>RECOMMENDED ANALYTICAL MASS</b>	<b>ESTIMATED IDL (<math>\mu\text{g/L}</math>)</b>
Aluminum	27	0.05
Antimony	121	0.08
Arsenic	75	0.9
Barium	137	0.5
Beryllium	9	0.1
Cadmium	111	0.1
Chromium	52	0.07
Cobalt	59	0.03
Copper	63	0.03
Lead	206, 207, 208	0.08
Manganese	55	0.1
Molybdenum	98	0.1
Nickel	60	0.2
Selenium	82	5
Silver	107	0.05
Thallium	205	0.09
Thorium	232	0.03
Uranium	238	0.02
Vanadium	51	0.02
Zinc	66	0.2

Instrument detection limits ( $3\sigma$ ) estimated from seven replicate integrations of the blank (1% v/v nitric acid) following calibration of the instrument with three replicate integrations of a multi-element standard.



TABLE 2: COMMON MOLECULAR ION INTERFERENCES IN ICP-MS

BACKGROUND MOLECULAR IONS

Molecular Ion	Mass	Element Interference <sup>a</sup>
NH <sup>+</sup>	15	
OH <sup>+</sup>	17	
OH <sub>2</sub> <sup>+</sup>	18	
C <sub>2</sub> <sup>+</sup>	24	
CN <sup>+</sup>	26	
CO <sup>+</sup>	28	
N <sub>2</sub> <sup>+</sup>	28	
N <sub>2</sub> H <sup>+</sup>	29	
NO <sup>+</sup>	30	
NOH <sup>+</sup>	31	
O <sub>2</sub> <sup>+</sup>	32	
O <sub>2</sub> H <sup>+</sup>	33	
<sup>36</sup> ArH <sup>+</sup>	37	
<sup>38</sup> ArH <sup>+</sup>	39	
<sup>40</sup> ArH <sup>+</sup>	41	
CO <sub>2</sub> <sup>+</sup>	44	
CO <sub>2</sub> H <sup>+</sup>	45	Sc
ArC <sup>+</sup> , ArO <sup>+</sup>	52	Cr
ArN <sup>+</sup>	54	Cr
ArNH <sup>+</sup>	55	Mn
ArO <sup>+</sup>	56	
ArOH <sup>+</sup>	57	
<sup>40</sup> Ar <sup>36</sup> Ar <sup>+</sup>	76	Se
<sup>40</sup> Ar <sup>38</sup> Ar <sup>+</sup>	78	Se
<sup>40</sup> Ar <sub>2</sub> <sup>+</sup>	80	Se

<sup>a</sup> method elements or internal standards affected by the molecular ions.

TABLE 2 (Continued).

MATRIX MOLECULAR IONS		
<b>CHLORIDE</b>		
Molecular Ion	Mass	Element Interference
$^{35}\text{ClO}^+$	51	V
$^{35}\text{ClOH}^+$	52	Cr
$^{37}\text{ClO}^+$	53	Cr
$^{37}\text{ClOH}^+$	54	Cr
$\text{Ar}^{35}\text{Cl}^+$	75	As
$\text{Ar}^{37}\text{Cl}^+$	77	Se
<b>SULPHATE</b>		
Molecular Ion	Mass	Element Interference
$^{32}\text{SO}^+$	48	
$^{32}\text{SOH}^+$	49	
$^{34}\text{SO}^+$	50	V, Cr
$^{34}\text{SOH}^+$	51	V
$\text{SO}_2^+, \text{S}_2^+$	64	Zn
$\text{Ar}^{32}\text{S}^+$	72	
$\text{Ar}^{34}\text{S}^+$	74	
<b>PHOSPHATE</b>		
Molecular Ion	Mass	Element Interference
$\text{PO}^+$	47	
$\text{POH}^+$	48	
$\text{PO}_2^+$	63	Cu
$\text{ArP}^+$	71	
<b>GROUP I, II METALS</b>		
Molecular Ion	Mass	Element Interference
$\text{ArNa}^+$	63	Cu
$\text{ArK}^+$	79	
$\text{ArCa}^+$	80	
<b>MATRIX OXIDES*</b>		
Molecular Ion	Masses	Element Interference
TiO	62-66	Ni, Cu, Zn
ZrO	106-112	Ag, Cd
MoO	108-116	Cd

\* Oxide interferences will normally be very small and will only impact the method elements when present at relatively high concentrations. Some examples of matrix oxides are listed of which the analyst should be aware. It is recommended that Ti and Zr isotopes are monitored in solid waste samples, which are likely to contain high levels of these elements. Mo is monitored as a method analyte.

**TABLE 3: INTERNAL STANDARDS AND LIMITATIONS OF USE**

Internal Standard	Mass	Possible Limitation
<b><sup>6</sup>Lithium</b>	6	<sup>a</sup>
<b>Scandium</b>	45	polyatomic ion interference
<b>Yttrium</b>	89	a,b
<b>Rhodium</b>	103	
<b>Indium</b>	115	isobaric interference by Sn
<b>Terbium</b>	159	
<b>Holmium</b>	165	
<b>Lutetium</b>	175	
<b>Bismuth</b>	209	a

a May be present in environmental samples.

b In some instruments Yttrium may form measurable amounts of  $YO^+$  (105 amu) and  $YOH^+$  (106 amu). If this is the case, care should be taken in the use of the cadmium elemental correction equation.

Internal standards recommended for use with this method are shown in bold face. Preparation procedures for these are included in section 7.3.

**TABLE 4: RECOMMENDED ANALYTICAL ISOTOPES AND ADDITIONAL  
MASSES WHICH MUST BE MONITORED**

Isotope	Element of Interest
<u>27</u>	Aluminum
<u>121, 123</u>	Antimony
<u>75</u>	Arsenic
<u>135, 137</u>	Barium
<u>9</u>	Beryllium
<u>106, 108, 111, 114</u>	Cadmium
<u>52, 53</u>	Chromium
<u>59</u>	Cobalt
<u>63, 65</u>	Copper
<u>206, 207, 208</u>	Lead
<u>55</u>	Manganese
<u>95, 97, 98</u>	Molybdenum
<u>60, 62</u>	Nickel
<u>77, 82</u>	Selenium
<u>107, 109</u>	Silver
<u>203, 205</u>	Thallium
<u>232</u>	Thorium
<u>238</u>	Uranium
<u>51</u>	Vanadium
<u>66, 67, 68</u>	Zinc
83	Krypton
99	Ruthenium
105	Palladium
118	Tin

NOTE: Isotopes recommended for analytical determination are underlined.

TABLE 5: RECOMMENDED ELEMENTAL EQUATIONS FOR DATA CALCULATIONS

Element	Elemental Equation	Note
Al	$(1.000)(^{27}\text{C})$	
Sb	$(1.000)(^{121}\text{C})$	
As	$(1.000)(^{75}\text{C}) - (3.127)[(^{77}\text{C}) - (0.815)(^{82}\text{C})]$	(1)
Ba	$(1.000)(^{137}\text{C})$	
Be	$(1.000)(^9\text{C})$	
Cd	$(1.000)(^{111}\text{C}) - (1.073)[(^{108}\text{C}) - (0.712)(^{106}\text{C})]$	(2)
Cr	$(1.000)(^{52}\text{C})$	(3)
Co	$(1.000)(^{59}\text{C})$	
Cu	$(1.000)(^{63}\text{C})$	
Pb	$(1.000)(^{206}\text{C}) + (1.000)(^{207}\text{C}) + (1.000)(^{208}\text{C})$	(4)
Mn	$(1.000)(^{55}\text{C})$	
Mo	$(1.000)(^{98}\text{C}) - (0.146)(^{99}\text{C})$	(5)
Ni	$(1.000)(^{60}\text{C})$	
Se	$(1.000)(^{82}\text{C})$	(6)
Ag	$(1.000)(^{107}\text{C})$	
Tl	$(1.000)(^{205}\text{C})$	
Th	$(1.000)(^{232}\text{C})$	
U	$(1.000)(^{238}\text{C})$	
V	$(1.000)(^{51}\text{C}) - (3.127)[(^{53}\text{C}) - (0.113)(^{52}\text{C})]$	(7)
Zn	$(1.000)(^{66}\text{C})$	

Cont.

TABLE 5 (Continued)

## INTERNAL STANDARDS

Element	Elemental Equation	Note
Bi	$(1.000)(^{209}\text{C})$	(8)
In	$(1.000)(^{115}\text{C}) - (0.016)(^{118}\text{C})$	
Sc	$(1.000)(^{45}\text{C})$	
Tb	$(1.000)(^{159}\text{C})$	
Y	$(1.000)(^{89}\text{C})$	

- 
- C - calibration blank subtracted counts at specified mass.
- (1) - correction for chloride interference with adjustment for Se77. ArCl 75/77 ratio may be determined from the reagent blank.
- (2) - correction for MoO interference. An additional isobaric elemental correction should be made if palladium is present.
- (3) - in 0.4% v/v HCl, the background from ClOH will normally be small. However the contribution may be estimated from the reagent blank.
- (4) - allowance for isotopic variability of lead isotopes.
- (5) - isobaric elemental correction for ruthenium.
- (6) - some argon supplies contain krypton as an impurity. Selenium is corrected for Kr82 by background subtraction.
- (7) - correction for chloride interference with adjustment for Cr53. ClO 51/53 ratio may be determined from the reagent blank.
- (8) - isobaric elemental correction for tin.

**TABLE 6: INSTRUMENT OPERATING CONDITIONS  
FOR PRECISION AND RECOVERY DATA**

---

Instrument	VG PlasmaQuad Type I
Plasma forward power	1.35 kW
Coolant flow rate	13.5 L/min
Auxiliary flow rate	0.6 L/min
Nebulizer flow rate	0.78 L/min
Solution uptake rate	0.6 mL/min
Spray chamber temperature	15°C

**Data Acquisition**

Detector mode	Pulse counting
Replicate integrations	3
Mass range	8 - 240 amu
Dwell time	320 $\mu$ s
Number of MCA channels	2048
Number of scan sweeps	85
Total acquisition time	3 minutes per sample

---

TABLE 7: TOTAL RECOVERABLE METHOD DETECTION LIMITS

ELEMENT	RECOMMENDED ANALYTICAL MASS	AQUEOUS μg/L	MDL* SOLIDS mg/kg
Aluminum	27	1.0	0.4
Antimony	121	0.4	0.2
Arsenic	75	1.4	0.6
Barium	137	0.8	0.4
Beryllium	9	0.3	0.1
Cadmium	111	0.5	0.2
Chromium	52	0.9	0.4
Cobalt	59	0.09	0.04
Copper	63	0.5	0.2
Lead	206, 207, 208	0.6	0.3
Manganese	55	0.1	0.05
Molybdenum	98	0.3	0.1
Nickel	60	0.5	0.2
Selenium	82	7.9	3.2
Silver	107	0.1	0.05
Thallium	205	0.3	0.1
Thorium	232	0.1	0.05
Uranium	238	0.1	0.05
Vanadium	51	2.5	1.0
Zinc	66	1.8	0.7

\* MDL concentrations are computed for original matrix with allowance for sample dilution during preparation.



**TABLE 8 : PRECISION AND RECOVERY DATA IN AQUEOUS MATRICES**  
**DRINKING WATER**

Element	Sample Concn. ( $\mu\text{g/L}$ )	Low Spike ( $\mu\text{g/L}$ )	Average Recovery R (%)	S(R)	RPD	High Spike ( $\mu\text{g/L}$ )	Average Recovery R (%)	S(R)	RPD
Al	175	50	115.8	5.9	0.4	200	102.7	1.6	1.1
Sb	<0.4	10	99.1	0.7	2.0	100	100.8	0.7	2.0
As	<1.4	50	99.7	0.8	2.2	200	102.5	1.1	2.9
Ba	43.8	50	94.8	3.9	5.8	200	95.6	0.8	1.7
Be	<0.3	10	113.5	0.4	0.9	100	111.0	0.7	1.8
Cd	<0.5	10	97.0	2.8	8.3	100	101.5	0.4	1.0
Cr	<0.9	10	111.0	3.5	9.0	100	99.5	0.1	0.2
Co	0.11	10	94.4	0.4	1.1	100	93.6	0.5	1.4
Cu	3.6	10	101.8	8.8	17.4	100	91.6	0.3	0.3
Pb	0.87	10	97.8	2.0	2.8	100	99.0	0.8	2.2
Mn	0.96	10	96.9	1.8	4.7	100	95.8	0.6	1.8
Mo	1.9	10	99.4	1.6	3.4	100	98.6	0.4	1.0
Ni	1.9	10	100.2	5.7	13.5	100	95.2	0.5	1.3
Se	<7.9	50	99.0	1.8	5.3	200	93.5	3.5	10.7
Ag	<0.1	50	100.7	1.5	4.2	200	99.0	0.4	1.0
Tl	<0.3	10	97.5	0.4	1.0	100	98.5	1.7	4.9
Th	<0.1	10	109.0	0.7	1.8	100	106.0	1.4	3.8
U	0.23	10	110.7	1.4	3.5	100	107.8	0.7	1.9
V	<2.5	50	101.4	0.1	0.4	200	97.5	0.7	2.1
Zn	5.2	50	103.4	3.3	7.7	200	96.4	0.5	1.0

S(R) Standard deviation of percent recovery.

RPD Relative percent difference between duplicate spike determinations.

< Sample concentration below established method detection limit.

TABLE 8 : PRECISION AND RECOVERY DATA IN AQUEOUS MATRICES (Cont).

WELL WATER

Element	Sample Conc. ( $\mu\text{g/L}$ )	Low Spike ( $\mu\text{g/L}$ )	Average Recovery R (%)	S(R)	RPD	High Spike ( $\mu\text{g/L}$ )	Average Recovery R (%)	S(R)	RPD
Al	34.3	50	100.1	3.9	0.8	200	102.6	1.1	1.3
Sb	0.46	10	98.4	0.9	1.9	100	102.5	0.7	1.9
As	<1.4	50	110.0	6.4	16.4	200	101.3	0.2	0.5
Ba	106	50	95.4	3.9	3.3	200	104.9	1.0	1.6
Be	<0.3	10	104.5	0.4	1.0	100	101.4	1.2	3.3
Cd	1.6	10	88.6	1.7	3.8	100	98.6	0.6	1.6
Cr	<0.9	10	111.0	0.0	0.0	100	103.5	0.4	1.0
Co	2.4	10	100.6	1.0	1.6	100	104.1	0.4	0.9
Cu	37.4	10	104.3	5.1	1.5	100	100.6	0.8	1.5
Pb	3.5	10	95.2	2.5	1.5	100	99.5	1.4	3.9
Mn	2770	10	*	*	1.8	100	*	*	0.7
Mo	2.1	10	103.8	1.1	1.6	100	102.9	0.7	1.9
Ni	11.4	10	116.5	6.3	6.5	100	99.6	0.3	0.0
Se	<7.9	50	127.3	8.4	18.7	200	101.3	0.2	0.5
Ag	<0.1	50	99.2	0.4	1.0	200	101.5	1.4	3.9
Tl	<0.3	10	93.9	0.1	0.0	100	100.4	1.8	5.0
Th	<0.1	10	103.0	0.7	1.9	100	104.5	1.8	4.8
U	1.8	10	106.0	1.1	1.6	100	109.7	2.5	6.3
V	<2.5	50	105.3	0.8	2.1	200	105.8	0.2	0.5
Zn	554	50	*	*	1.2	200	102.1	5.5	3.2

S(R) Standard deviation of percent recovery.  
 RPD Relative percent difference between duplicate spike determinations.  
 < Sample concentration below established method detection limit.  
 \* Spike concentration <10% of sample background concentration.

TABLE 8 : PRECISION AND RECOVERY DATA IN AQUEOUS MATRICES (Cont).

POND WATER

Element	Sample Concn. (µg/L)	Low Spike (µg/L)	Average Recovery R (%)	S(R)	RPD	High Spike (µg/L)	Average Recovery R (%)	S(R)	RPD
Al	610	50	*	*	1.7	200	78.2	9.2	5.5
Sb	<0.4	10	101.1	1.1	2.9	100	101.5	3.0	8.4
As	<1.4	50	100.8	2.0	5.6	200	96.8	0.9	2.6
Ba	28.7	50	102.1	1.8	2.4	200	102.9	3.7	9.0
Be	<0.3	10	109.1	0.4	0.9	100	114.4	3.9	9.6
Cd	<0.5	10	106.6	3.2	8.3	100	105.8	2.8	7.6
Cr	2.0	10	107.0	1.0	1.6	100	100.0	1.4	3.9
Co	0.79	10	101.6	1.1	2.7	100	101.7	1.8	4.9
Cu	5.4	10	107.5	1.4	1.9	100	98.1	2.5	6.8
Pb	1.9	10	108.4	1.5	3.2	100	106.1	0.0	0.0
Mn	617	10	*	*	1.1	100	139.0	11.1	4.0
Mo	0.98	10	104.2	1.4	3.5	100	104.0	2.1	5.7
Ni	2.5	10	102.0	2.3	4.7	100	102.5	2.1	5.7
Se	<7.9	50	102.7	5.6	15.4	200	105.5	1.4	3.8
Ag	0.12	50	102.5	0.8	2.1	200	105.2	2.7	7.1
Tl	<0.3	10	108.5	3.2	8.3	100	105.0	2.8	7.6
Th	0.19	10	93.1	3.5	10.5	100	93.9	1.6	4.8
U	0.30	10	107.0	2.8	7.3	100	107.2	1.8	4.7
V	3.5	50	96.1	5.2	14.2	200	101.5	0.2	0.5
Zn	6.8	50	99.8	1.7	3.7	200	100.1	2.8	7.7

S(R) Standard deviation of percent recovery.  
 RPD Relative percent difference between duplicate spike determinations.  
 < Sample concentration below established method detection limit.  
 \* Spike concentration <10% of sample background concentration.

TABLE 8 : PRECISION AND RECOVERY DATA IN AQUEOUS MATRICES (Cont).

SEWAGE TREATMENT PRIMARY EFFLUENT

Element	Sample Concn. ( $\mu\text{g/L}$ )	Low Spike ( $\mu\text{g/L}$ )	Average Recovery R (%)	S(R)	RPD	High Spike ( $\mu\text{g/L}$ )	Average Recovery R (%)	S(R)	RPD
Al	1150	50	*	*	3.5	200	100.0	13.8	1.5
Sb	1.5	10	95.7	0.4	0.9	100	104.5	0.7	1.9
As	<1.4	50	104.2	4.5	12.3	200	101.5	0.7	2.0
Ba	202	50	79.2	9.9	2.5	200	108.6	4.6	5.5
Be	<0.3	10	110.5	1.8	4.5	100	106.4	0.4	0.9
Cd	9.2	10	101.2	1.3	0.0	100	102.3	0.4	0.9
Cr	128	10	*	*	1.5	100	102.1	1.7	0.4
Co	13.4	10	95.1	2.7	2.2	100	99.1	1.1	2.7
Cu	171	10	*	*	2.4	100	105.2	7.1	0.7
Pb	17.8	10	95.7	3.8	1.1	100	102.7	1.1	2.5
Mn	199	10	*	*	1.5	100	103.4	2.1	0.7
Mo	136	10	*	*	1.4	100	105.7	2.4	2.1
Ni	84.0	10	88.4	16.3	4.1	100	98.0	0.9	0.0
Se	<7.9	50	112.0	10.9	27.5	200	108.8	3.0	7.8
Ag	10.9	50	97.1	0.7	1.5	200	102.6	1.4	3.7
Tl	<0.3	10	97.5	0.4	1.0	100	102.0	0.0	0.0
Th	0.11	10	15.4	1.8	30.3	100	29.3	0.8	8.2
U	0.71	10	109.4	1.8	4.3	100	109.3	0.7	1.8
V	<2.5	50	90.9	0.9	0.6	200	99.4	2.1	6.0
Zn	163	50	85.8	3.3	0.5	200	102.0	1.5	1.9

S(R) Standard deviation of percent recovery.  
 RPD Relative percent difference between duplicate spike determinations.  
 < Sample concentration below established method detection limit.  
 \* Spike concentration <10% of sample background concentration.

TABLE 8 : PRECISION AND RECOVERY DATA IN AQUEOUS MATRICES (Cont).

INDUSTRIAL EFFLUENT

Element	Sample Concn. ( $\mu\text{g/L}$ )	Low Spike ( $\mu\text{g/L}$ )	Average Recovery R (%)	S(R)	RPD	High Spike ( $\mu\text{g/L}$ )	Average Recovery R (%)	S(R)	RPD
Al	44.7	50	98.8	8.7	5.7	200	90.4	2.1	2.2
Sb	2990	10	*	*	0.3	100	*	*	0.0
As	<1.4	50	75.1	1.8	6.7	200	75.0	0.0	0.0
Ba	100	50	96.7	5.5	3.4	200	102.9	1.1	0.7
Be	<0.3	10	103.5	1.8	4.8	100	100.0	0.0	0.0
Cd	10.1	10	106.5	4.4	2.4	100	97.4	1.1	2.8
Cr	171	10	*	*	0.0	100	127.7	2.4	1.7
Co	1.3	10	90.5	3.2	8.7	100	90.5	0.4	1.3
Cu	101	10	*	*	0.9	100	92.5	2.0	1.6
Pb	294	10	*	*	2.6	100	108.4	2.1	0.0
Mn	154	10	*	*	2.8	100	103.6	3.7	1.6
Mo	1370	10	*	*	1.4	100	*	*	0.7
Ni	17.3	10	107.4	7.4	5.0	100	88.2	0.7	1.0
Se	15.0	50	129.5	9.3	15.1	200	118.3	1.9	3.6
Ag	<0.1	50	91.8	0.6	1.7	200	87.0	4.9	16.1
Tl	<0.3	10	90.5	1.8	5.5	100	98.3	1.0	2.8
Th	0.29	10	109.6	1.2	2.7	100	108.7	0.0	0.0
U	0.17	10	104.8	2.5	6.6	100	109.3	0.4	0.9
V	<2.5	50	74.9	0.1	0.3	200	72.0	0.0	0.0
Zn	43.4	50	85.0	4.0	0.6	200	97.6	1.0	0.4

S(R) Standard deviation of percent recovery.

RPD Relative percent difference between duplicate spike determinations.

< Sample concentration below established method detection limit.

\* Spike concentration <10% of sample background concentration.

TABLE 9 : PRECISION AND RECOVERY DATA IN SOLID MATRICES

EPA HAZARDOUS SOIL #884

Element	Sample Concn. (mg/kg)	Low+ Spike (mg/kg)	Average Recovery R (%)	S(R)	RPD	High+ Spike (mg/kg)	Average Recovery R (%)	S(R)	RPD
Al	5170	20	*	*	-	100	*	*	-
Sb	5.4	20	69.8	2.5	4.7	100	70.4	1.8	6.5
As	8.8	20	104.7	5.4	9.1	100	102.2	2.2	5.4
Ba	113	20	54.9	63.6	18.6	100	91.0	9.8	0.5
Be	0.6	20	100.1	0.6	1.5	100	102.9	0.4	1.0
Cd	1.8	20	97.3	1.0	1.4	100	101.7	0.4	1.0
Cr	83.5	20	86.7	16.1	8.3	100	105.5	1.3	0.0
Co	7.1	20	98.8	1.2	1.9	100	102.9	0.7	1.8
Cu	115	20	86.3	13.8	3.4	100	102.5	4.2	4.6
Pb	152	20	85.0	45.0	13.9	100	151.7	25.7	23.7
Mn	370	20	*	*	12.7	100	85.2	10.4	2.2
Mo	4.8	20	95.4	1.5	2.9	100	95.2	0.7	2.0
Ni	19.2	20	101.7	3.8	1.0	100	102.3	0.8	0.8
Se	<3.2	20	79.5	7.4	26.4	100	100.7	9.4	26.5
Ag	1.1	20	96.1	0.6	0.5	100	94.8	0.8	2.3
Tl	0.24	20	94.3	1.1	3.1	100	97.9	1.0	2.9
Th	1.0	20	69.8	0.6	1.3	100	76.0	2.2	7.9
U	1.1	20	100.1	0.2	0.0	100	102.9	0.0	0.0
V	17.8	20	109.2	4.2	2.3	100	106.7	1.3	2.4
Zn	128	20	87.0	27.7	5.5	100	113.4	12.9	14.1

S(R) Standard deviation of percent recovery.  
 RPD Relative percent difference between duplicate spike determinations.  
 < Sample concentration below established method detection limit.  
 \* Spike concentration <10% of sample background concentration.  
 - Not determined.  
 + Equivalent.

TABLE 9 : PRECISION AND RECOVERY DATA IN SOLID MATRICES (Cont).

NBS 1645 RIVER SEDIMENT

Element	Sample Concn. (mg/kg)	Low+ Spike (mg/kg)	Average Recovery R (%)	S(R)	RPD	High+ Spike (mg/kg)	Average Recovery R (%)	S(R)	RPD
Al	5060	20	*	*	-	100	*	*	-
Sb	21.8	20	73.9	6.5	9.3	100	81.2	1.5	3.9
As	67.2	20	104.3	13.0	7.6	100	107.3	2.1	2.9
Ba	54.4	20	105.6	4.9	2.8	100	98.6	2.2	3.9
Be	0.59	20	88.8	0.2	0.5	100	87.9	0.1	0.2
Cd	8.3	20	92.9	0.4	0.0	100	95.7	1.4	3.9
Cr	29100	20	*	*	-	100	*	*	-
Co	7.9	20	97.6	1.3	2.6	100	103.1	0.0	0.0
Cu	112	20	121.0	9.1	1.5	100	105.2	2.2	1.8
Pb	742	20	*	*	-	100	-	-	-
Mn	717	20	*	*	-	100	-	-	-
Mo	17.1	20	89.8	8.1	12.0	100	98.4	0.7	0.9
Ni	41.8	20	103.7	6.5	4.8	100	102.2	0.8	0.0
Se	<3.2	20	108.3	14.3	37.4	100	93.9	5.0	15.1
Ag	1.8	20	94.8	1.6	4.3	100	96.2	0.7	1.9
Tl	1.2	20	91.2	1.3	3.6	100	94.4	0.4	1.3
Th	0.90	20	91.3	0.9	2.6	100	92.3	0.9	2.8
U	0.79	20	95.6	1.8	5.0	100	98.5	1.2	3.5
V	21.8	20	91.8	4.6	5.7	100	100.7	0.6	0.8
Zn	1780	20	*	*	-	100	*	*	-

S(R) Standard deviation of percent recovery.  
 RPD Relative percent difference between duplicate spike determinations.  
 < Sample concentration below established method detection limit.  
 \* Spike concentration <10% of sample background concentration.  
 - Not determined.  
 + Equivalent.

TABLE 9 : PRECISION AND RECOVERY DATA IN SOLID MATRICES (Cont).

EPA ELECTROPLATING SLUDGE #286

Element	Sample Concn. (mg/kg)	Low+ Spike (mg/kg)	Average Recovery R (%)	S(R)	RPD	High+ Spike (mg/kg)	Average Recovery R (%)	S(R)	RPD
Al	5110	20	*	*	-	100	*	*	-
Sb	8.4	20	55.4	1.5	4.1	100	61.0	0.2	0.9
As	41.8	20	91.0	2.3	1.7	100	94.2	0.8	1.5
Ba	27.3	20	1.8	7.1	8.3	100	0	1.5	10.0
Be	0.25	20	92.0	0.9	2.7	100	93.4	0.3	0.9
Cd	112	20	85.0	5.2	1.6	100	88.5	0.8	0.5
Cr	7980	20	*	*	-	100	*	*	-
Co	4.1	20	89.2	1.8	4.6	100	88.7	1.5	4.6
Cu	740	20	*	*	6.0	100	61.7	20.4	5.4
Pb	1480	20	*	*	-	100	*	*	-
Mn	295	20	*	*	-	100	-	-	-
Mo	13.3	20	82.9	1.2	1.3	100	89.2	0.4	1.0
Ni	450	20	*	*	6.8	100	83.0	10.0	4.5
Se	3.5	20	89.7	3.7	4.2	100	91.0	6.0	18.0
Ag	5.9	20	89.8	2.1	4.6	100	85.1	0.4	1.1
Tl	1.9	20	96.9	0.9	2.4	100	98.9	0.9	2.4
Th	3.6	20	91.5	1.3	3.2	100	97.4	0.7	2.0
U	2.4	20	107.7	2.0	4.6	100	109.6	0.7	1.8
V	21.1	20	105.6	1.8	2.1	100	97.4	1.1	2.5
Zn	13300	20	*	*	-	100	*	*	-

S(R) Standard deviation of percent recovery.  
 RPD Relative percent difference between duplicate spike determinations.  
 < Sample concentration below established method detection limit.  
 \* Spike concentration <10% of sample background concentration.  
 - Not determined.  
 + Equivalent.